


Upgrade and Migration Guide

Oracle[®] Health Sciences InForm 6.0.0.2

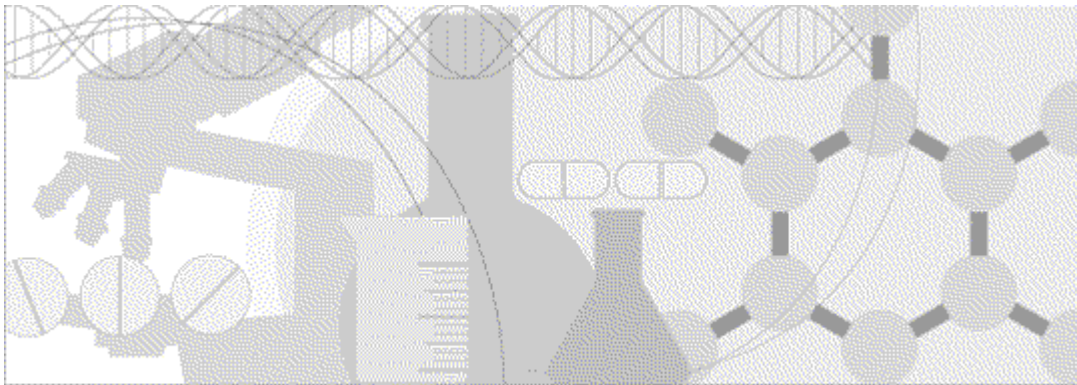

ORACLE[®]

Copyright © 2014, Oracle and/or its affiliates. All rights reserved.

The Programs (which include both the software and documentation) contain proprietary information; they are provided under a license agreement containing restrictions on use and disclosure and are also protected by copyright, patent, and other intellectual and industrial property laws. Reverse engineering, disassembly, or decompilation of the Programs, except to the extent required to obtain interoperability with other independently created software or as specified by law, is prohibited.

The information contained in this document is subject to change without notice. If you find any problems in the documentation, please report them to us in writing. This document is not warranted to be error-free. Except as may be expressly permitted in your license agreement for these Programs, no part of these Programs may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose.

If the Programs are delivered to the United States Government or anyone licensing or using the Programs on behalf of the United States Government, the following notice is applicable:

U.S. GOVERNMENT RIGHTS Programs, software, databases, and related documentation and technical data delivered to U.S. Government customers are "commercial computer software" or "commercial technical data" pursuant to the applicable Federal Acquisition Regulation and agency-specific supplemental regulations. As such, use, duplication, disclosure, modification, and adaptation of the Programs, including documentation and technical data, shall be subject to the licensing restrictions set forth in the applicable Oracle license agreement, and, to the extent applicable, the additional rights set forth in FAR 52.227-19, Commercial Computer Software -- Restricted Rights (June 1987). Oracle USA, Inc., 500 Oracle Parkway, Redwood City, CA 94065.

The Programs are not intended for use in any nuclear, aviation, mass transit, medical, or other inherently dangerous applications. It shall be the licensee's responsibility to take all appropriate fail-safe, backup, redundancy and other measures to ensure the safe use of such applications if the Programs are used for such purposes, and we disclaim liability for any damages caused by such use of the Programs.

The Programs may provide links to Web sites and access to content, products, and services from third parties. Oracle is not responsible for the availability of, or any content provided on, third-party Web sites. You bear all risks associated with the use of such content. If you choose to purchase any products or services from a third party, the relationship is directly between you and the third party. Oracle is not responsible for: (a) the quality of third-party products or services; or (b) fulfilling any of the terms of the agreement with the third party, including delivery of products or services and warranty obligations related to purchased products or services. Oracle is not responsible for any loss or damage of any sort that you may incur from dealing with any third party.

Oracle is a registered trademark of Oracle Corporation and/or its affiliates. Other names may be trademarks of their respective owners.

This documentation may include references to materials, offerings, or products that were previously offered by Phase Forward Inc. Certain materials, offerings, services, or products may no longer be offered or provided. Oracle and its affiliates cannot be held responsible for any such references should they appear in the text provided.

Contents

About this guide vii

Overview of this guide.....	viii
Audience	viii
Documentation	ix
If you need assistance.....	xi

Chapter 1 Upgrading and migrating to the InForm 6.0.0.2 release 1

Overview of upgrading and migrating to this release.....	2
Requirements for the InForm studies you can migrate	3
Hardware and software requirements for the InForm 6.0.0.2 release.....	4

Chapter 2 Upgrading the database environment 5

Overview of upgrading the database environment.....	6
---	---

Chapter 3 Upgrading and migrating the InForm studies to this release 7

Overview of upgrading and migrating InForm studies to this release.....	8
Upgrade to a new Oracle database environment	9
Preparing to upgrade or migrate the InForm software and studies.....	10
Considerations before you migrate	11
Changed resource files.....	15
Create backup files	19
Prepare the InForm 4.6 study database schema for backup.....	20
Back up the InForm study database.....	21
Record study, customer-defined database, and connection information	23
Create a Central Designer 2.0 deployment package	23
Migrating InForm studies to a new machine	24
Install the InForm 6.0.0.2 software.....	25
Create the database user.....	26
Create tablespaces other than INFORM	26
Import the DMP files	27
Grant the database user privileges.....	29
Set up the InForm servers, studies, and customer-defined databases	29
Upgrade the study schemas to the InForm 6.0.0.2 release	30
Update the statistics for the InForm studies	31
Start the InForm server.....	31
Create new database objects for the review schema	31
Set the MEDMLINSTALLER server.....	32
Restore the Central Designer rule engines to the original location	32
Run the 64-bit Central Designer deployment package	33
Upgrade the study resources	33
Terminate obsolete users	34
Set the locale for InForm 4.6 sites and users.....	35
Delete PF_PATH table records.....	36
Perform the final study configuration.....	36
Restore the Windows registry settings.....	37
Install or upgrade the InForm Adapter software.....	37
Start the InForm studies.....	37
Running an in-place upgrade of any InForm 6.0 release	38
Post-upgrade tasks	40

Clear the client-side browser caches	40
Verify the rights assigned to each rights group	41
Reset your password	42

Chapter 4 Upgrading and migrating the reporting environment to this release 43

Overview of upgrading and migrating the reporting environment to this release	44
Migrating the reporting environment to a new machine	45
Install and configure the Cognos 8.4.1 environment	46
Create an export package and export the data from your Cognos environment	48
Create the streams between the InForm study schema and the Reporting schema, and install the Reporting schema	50
Import the export package	51
Run the PFCognosConfig utility	55
Modify the PFRSetupTrial.xml file	55
Run PFRInit on the InForm application server	56
Delete obsolete reports	56
Associate custom reports with the new operational packages	57
Restart the InForm Model Updater service	57
Revalidate custom reports	58
Update schedules and jobs	58
Delete the study organizational units from the LDAP configuration	58
Running an in-place upgrade of an InForm 6.0 reporting environment	59

Chapter 5 Migrating the InForm Portal software 61

Prepare the InForm Portal 4.6 database schema for backup	62
Back up the InForm Portal database and customized files	63
Back up the InForm Portal 4.6 or 5.0 database and customized files—Oracle 10g	63
Back up the InForm Portal 4.6, 5.5, or 6.0 database and customized files—Oracle 11g	64
Migrate the InForm Portal software to a new database instance	65
Migrate the InForm Portal software in the same database instance	68
Reset the InForm Portal Home page	69

Chapter 6 Database and file changes 71

Database schema changes	72
Changes to the database schema from the InForm 4.6 software to the InForm 4.7 software	72
Changes to the database schema from the InForm 4.7 software to the InForm 5.0 software	73
Changes to the database schema from the InForm 5.0 software to the InForm 5.5 software	74
Changes to the database schema from the InForm 5.5 software to the InForm 6.0 software	75
Data changes	79
Data changes from the InForm 4.6 software to the InForm 4.7 software	79
Data changes from the InForm 4.7 software to the InForm 5.0 software	79
Data changes from the InForm 5.0 software to the InForm 5.5 software	80
Data changes from the InForm 5.0 software to the InForm 5.5 software	80
MedML schema changes	81
MedML schema changes from the InForm 4.6 software to the InForm 4.7 software	81
MedML schema changes from the InForm 4.7 software to the InForm 5.0 software	82
MedML schema changes from the InForm 5.0 software to the InForm 5.5 software	82
MedML schema changes from the InForm 5.5 software to the InForm 6.0 software	83
Base MedML changes	85
Base MedML changes from the InForm 4.6 software to the InForm 4.7 software	85
Base MedML changes from the InForm 4.7 software to the InForm 5.0 software	87
Base MedML changes from the InForm 5.0 software to the InForm 5.5 software	88
Base MedML changes from the InForm 5.5 software to the InForm 6.0 software	89
Design changes	91

Design changes for the InForm 6.0 software.....	91
Database and file changes to support the InForm Portal migration	92
Database schema changes for the InForm Portal software.....	92
Data changes for the InForm Portal.....	92
Additions to the PF_CATEGORIES table for the InForm Portal software.....	92
Records inserted into the PF_SUBCAT table for the InForm Portal software.....	93
File changes to support the InForm Portal software migration	94

About this guide

In this preface

Overview of this guideviii

Documentationix

If you need assistance.....xi

Overview of this guide

The *Upgrade and Migration Guide* provides instructions for upgrading and migrating the InForm software and InForm Portal software to the current InForm release, and for upgrading the Cognos 8 Business Intelligence software for use with the Reporting and Analysis module. The guide also describes any changes and additions made to the database schema, MedML, and resource files.

Audience

This guide is for database and system administrators who are responsible for upgrading or migrating the InForm software and InForm Portal software to the current InForm release, and for upgrading the Cognos 8 Business Intelligence software for use with the Reporting and Analysis module.

Documentation

All documentation is available from the Oracle Software Delivery Cloud (<https://edelivery.oracle.com>) and the Download Center (<https://extranet.phaseforward.com>).

Document	Description
<i>Release Notes</i>	The <i>Release Notes</i> document describes enhancements introduced and problems fixed in the current release, upgrade considerations, release history, and other late-breaking information.
<i>Known Issues</i>	The <i>Known Issues</i> document provides detailed information about the known issues in this release, along with workarounds, if available.
<i>Upgrade and Migration Guide</i>	The <i>Upgrade and Migration Guide</i> provides instructions for upgrading and migrating the InForm software and InForm Portal software to the current InForm release, and for upgrading the Cognos 8 Business Intelligence software for use with the Reporting and Analysis module. The guide also describes any changes and additions made to the database schema, MedML, and resource files.
<i>Secure Configuration Guide</i>	The <i>Secure Configuration Guide</i> provides an overview of the security features provided with the Oracle® Health Sciences InForm application, including details about the general principles of application security, and how to install, configure, and use the InForm application securely.
<i>Installation Guide</i>	The <i>Installation Guide</i> describes how to install the software and configure the environment for the InForm application and Cognos 8 Business Intelligence software.
<i>Study and Reporting Setup Guide</i>	The <i>Study and Reporting Setup Guide</i> describes how to perform the tasks that are required to set up an InForm study and configure the Reporting and Analysis module for the study.
<i>User Guide</i>	<p>The <i>User Guide</i> provides an overview of the InForm application including details on multilingual studies, how to navigate through the user interface, and how to use the application to accomplish typical tasks you perform while running a clinical study.</p> <p>This document is also available from the Documentation CD and the InForm user interface.</p>
<i>Reporting and Analysis Guide</i>	<p>The <i>Reporting and Analysis Guide</i> provides an overview of the Reporting and Analysis module. It includes a brief overview of the Reporting and Analysis interface, illustrates how to access the Ad Hoc Reporting feature, and describes the study management and clinical data packages available for Reporting and Analysis. It also provides detailed descriptions of each standard report that is included with your installation.</p> <p>This document is also available from the Documentation CD.</p>
<i>Reporting Database Schema Guide</i>	The <i>Reporting Database Schema Guide</i> describes the Reporting and Analysis database schema, and provides information on creating Reporting Database Extracts (RDEs).

Document	Description
<i>Portal Administration Guide</i>	<p>The <i>Portal Administration Guide</i> provides step-by-step instructions for setting up the InForm Portal software, and configuring and managing the InForm Portal application.</p> <p>This document is also available from the Documentation CD.</p>
<i>Utilities Guide</i>	<p>The <i>Utilities Guide</i> provides information about and step-by-step instructions for using the following utilities:</p> <ul style="list-style-type: none"> • PFConsole utility • MedML Installer utility • InForm Data Import utility • InForm Data Export utility • InForm Performance Monitor utility • InForm Report Folder Maintenance utility <p>This guide also provides reference information for the MedML elements and scripting objects that are used to import and export data to and from the InForm application, as well as sample data import XML.</p> <p>This document is also available from the Documentation CD.</p>
MedML Installer utility online Help	<p>The MedML Installer utility online Help provides information about, and step-by-step instructions for using, the MedML Installer utility, which is used to load XML that defines study components into the InForm database.</p> <p>This guide also provides reference information for the MedML elements and scripting objects that are used to import and export data to and from the InForm application, as well as sample data import XML.</p> <p>This document is also available from the user interface.</p>
InForm Data Export utility online Help	<p>The InForm Data Export utility online Help provides information about and step-by-step instructions for using the InForm Data Export utility, which is used to export data from the InForm application to the following output formats:</p> <ul style="list-style-type: none"> • Customer-defined database (CDD). • Name value pairs. <p>This document is also available from the user interface.</p>
InForm Data Import utility online Help	<p>The InForm Data Import utility online Help provides information about and step-by-step instructions for using the InForm Data Import utility, which is used to import data into the InForm application.</p> <p>This document is also available from the user interface.</p>

If you need assistance

Oracle customers have access to support through My Oracle Support. For information, visit <http://www.oracle.com/pls/topic/lookup?ctx=acc&id=info>, or if you are hearing impaired, visit <http://www.oracle.com/pls/topic/lookup?ctx=acc&id=trs>.

CHAPTER 1

Upgrading and migrating to the InForm 6.0.0.2 release

In this chapter

Overview of upgrading and migrating to this release	2
Requirements for the InForm studies you can migrate.....	3
Hardware and software requirements for the InForm 6.0.0.2 release.....	4

Overview of upgrading and migrating to this release

You can perform a new installation of the InForm 6.0.0.2 software, or you can upgrade and migrate to the InForm 6.0.0.2 software from a previous version of the InForm software.

To perform a new installation, see the *Installation Guide*.

To perform an upgrade and migration, perform the following steps, as needed:

- **Upgrade the database environment**—If you are upgrading to the InForm 6.0.0.2 release from an environment that uses the Oracle 10g database or Oracle 11g database software. For more information, see *Upgrading the database environment* (on page 5).
- **Migrate the InForm studies**—You can upgrade to the InForm 6.0.0.2 software from the InForm 4.6 SP2 release or later, any InForm 5.0, 5.5, or 6.0 release.

- You can migrate any InForm 4.6, 5.0, 5.5, or 6.0 study to a new machine using a series of scripts.

Oracle does not support an in-place upgrade for any InForm 4.6, 5.0, or 5.5 release.

- You can migrate any InForm 6.0 release to a new machine, or you can perform an in-place upgrade.

For more information, see *Migrating InForm studies to a new machine* (on page 24).

- **Upgrade the reporting environment**—Optional. Use this step to upgrade your reporting environment for any InForm 4.6 SP2 studies or later, any InForm 5.0, 5.5 studies, or 6.0 studies. This is required only if you are using the Reporting and Analysis module.

- You can migrate the reporting environment for any InForm 4.6, 5.0, 5.5, or 6.0 study to a new machine.

Oracle does not support an in-place upgrade for any InForm 4.6, InForm 5.0, or InForm 5.5 release. You must install a new Cognos 8.4.1 64-bit environment, and create, export, and then import your export packages.

- For InForm 6.0 releases, you can migrate to a new machine or you can perform an in-place upgrade.

For more information, see *Upgrading and migrating the reporting environment to this release* (on page 43).

- **Migrate the InForm Portal software**—Optional. Use this step to upgrade or migrate the InForm Portal software for the InForm 6.0.0.2 release. This is required only if you are using the InForm Portal software. For more information, see *Migrating the InForm Portal software* (on page 61).

Requirements for the InForm studies you can migrate

You can upgrade to the InForm 6.0.0.2 release from the InForm 4.6 SP2 release or later, any InForm 5.0, 5.5, or 6.0 releases.

For any InForm 4.6 SP2 release or later, or any InForm 5.0 or 5.5 study, you must create a new deployment package with a new study version in the Central Designer 2.0 application. For more information, see ***Create a Central Designer 2.0 deployment package*** (on page 23). This step is not required if you are upgrading from any InForm 6.0 release.

There are no limitations on upgrading or migrating studies from any InForm 5.0, 5.5, or 6.0 release.

However, you can only migrate InForm 4.6 studies that meet the following criteria.

- The study is deployed in the InForm 4.6 SP2 release or later.
- The study must be deployed in the English (en-US) language.

The Study Version in the Central Designer application indicates the study language.

You cannot migrate studies to the InForm 6.0.0.2 release that:

- Use synchronization.
- Were deployed in any of the following InForm releases, even if they were updated to the InForm 4.6 SP2 release:
 - InForm 4.0
 - InForm 4.1
 - InForm 4.5 prior to the InForm 4.5 SP3b release.

Hardware and software requirements for the InForm 6.0.0.2 release

Before upgrading and migrating to the InForm 6.0.0.2 release, review the hardware and software requirements in the *Release Notes*.

CHAPTER 2

Upgrading the database environment

In this chapter

Overview of upgrading the database environment	6
--	---

Overview of upgrading the database environment

If you are upgrading to the InForm 6.0.0.2 release from an environment that uses the Oracle 10g database software, you must upgrade to the Oracle 11g database software.

- For information on which versions of the Oracle database software are supported, see the *Release Notes*.
- For information on installing the Oracle 11g database software, see your Oracle database documentation.
- For information on the database parameter settings, and guidance on installing and configuring the Oracle database software for the InForm environment, see the *Installation Guide*.

Note: Oracle does not support an in-place or direct upgrade from the Oracle 10g database to the Oracle 11g database software. You must export studies from an Oracle 10g environment, import them into an Oracle 11g environment, and then proceed with the rest of the upgrade.

CHAPTER 3

Upgrading and migrating the InForm studies to this release

In this chapter

Overview of upgrading and migrating InForm studies to this release..... 8

Upgrade to a new Oracle database environment..... 9

Preparing to upgrade or migrate the InForm software and studies..... 10

Migrating InForm studies to a new machine..... 24

Running an in-place upgrade of any InForm 6.0 release..... 38

Post-upgrade tasks 40

Overview of upgrading and migrating InForm studies to this release

You can upgrade to the InForm 6.0.0.2 software from the InForm 4.6 SP2 release or later, or any InForm 5.0, 5.5, or 6.0 releases using a series of script.

For more information, see:

- ***Requirements for the InForm studies you can migrate*** (on page 3).
- ***Migrating InForm studies to a new machine*** (on page 24).

For InForm 6.0 releases, you can migrate to a new machine or you can perform an in-place upgrade.

For more information, see ***Running an in-place upgrade of any InForm 6.0 release*** (on page 38). Oracle does not support an in-place upgrade for any InForm 4.6, 5.0, or 5.5 release.

Upgrade to a new Oracle database environment

- 1 Create a new Oracle 11.2 database environment.
- 2 Create new study databases in the new Oracle 11.2 environment.

For more information, see your Oracle database documentation.

Preparing to upgrade or migrate the InForm software and studies

To prepare to migrate the InForm studies:

- 1 Review the following topics:
 - *Database and file customizations* (on page 11).
 - *Character set differences* (on page 12).
 - *Query length limitations* (on page 12).
 - *Oracle 11g database import recommendations* (on page 13).
 - *Rights assigned to right groups* (on page 13).
 - *Itemset types and migrated studies* (on page 14).
 - *Imported child controls and migrated studies* (on page 14).
 - *Changed resource files* (on page 15).
- 2 *Create backup files* (on page 19).
- 3 *Prepare the InForm 4.6 study database schema for backup* (on page 20).
- 4 Back up the InForm study database:
 - *Back up the InForm 4.6 or 5.0 database and customized files—Oracle 10g* (on page 21).
 - *Back up the InForm 4.6, 5.5, or 6.0 database and customized files—Oracle 11g* (on page 22).
- 5 *Record study, customer-defined database, and connection information* (on page 23).
- 6 *Create a Central Designer 2.0 deployment package* (on page 23).

Note: This step is not required if you are upgrading from any InForm 6.0 release.

Considerations before you migrate

Database and file customizations

For lists of changes to the InForm 6.0.0.2 software that might affect your customizations, see *Database and file changes* (on page 71).

- Ensure that every tablespace in the study schema that includes InForm tables exists on the target database server before you import the study schema.
- Database schema objects created or customized outside of the standard InForm application are not supported for upgrade and can cause errors during the migration process.

Note: Oracle does not support migrating tables created outside of the standard InForm application schema or by Enterprise Adoption customers. These tables might cause errors when you import the study DMP files.

- Any schema customization must have an upgrade path that results in a valid InForm 6.0 schema at the point in the migration where you create the review schema. For more information, see *Create new database objects for the review schema* (on page 31).

You can record the customizations, and reintroduce them after you complete the standard migration procedures. Such customizations might include changes to:

- Views
- Functions
- Procedures
- Packages
- Triggers

Note: If any schema objects such as views, functions, procedures, packages, or triggers are invalid, pfadmin RECREATEREVIEWSCHEMA fails. You must repair any invalid schema objects before you can rerun pfadmin RECREATEREVIEWSCHEMA.

- Any study resources that have been customized for deployment must be addressed separately before migration.

The customized changes must be modified to comply with the InForm 6.0 software and InForm 6.0 MedML schema and reinstalled so that the customized changes are not lost due to the upgrade.

Character set differences

Character set differences between the InForm 4.6 release and the InForm 6.0.0.2 release might cause issues with table column widths.

To prevent these issues, you must prepare the InForm 4.6 database schemas before you migrate.

For more information, see:

- ***Prepare the InForm 4.6 study database schemas for backup*** (on page 20).
- ***Prepare the InForm 4.6 InForm Portal database schema for backup*** (on page 62).

Query length limitations

If you are migrating a study from an InForm 4.6 release, you need to review the query lengths.

If a query includes Latin-1 characters and is close to the query text size limit (4000 bytes), the query may exceed the size limit after the Latin-1 characters are converted to UTF-8 characters during the migration to the InForm 6.0.0.2 release.

If the query size exceeds the field limit, the query will be truncated.

Contact Oracle Global Support if you find any queries that have been truncated during the migration.

Oracle 11g database import recommendations

Oracle recommends the following when importing previously exported studies into an Oracle 11g database.

- Use single-user imports for importing study data into the new Oracle database.

Note: Do not import streams metadata from an Oracle 10g database to an Oracle 11g database.

Before you import the study into the study database, set the Oracle database import parameters **STREAMS_CONFIGURATION** and **STREAMS_INSTANTIATION** to **N**.

- Do not import the **STRMADMIN** user and **STRMADMIN** proxy user into the Oracle 11g database.
- Set **grants=N** for the import.

If you set **grants=Y**, grants to the streams proxy user from the Oracle 10g database may be attempted, but might not work, causing error messages to appear during the import.

Note: If you are importing study databases that were used in an environment that includes the InForm Adapter application, you might encounter **Oracle error 31167: XML nodes over 64K in size cannot be inserted**. You can ignore this error.

Note: If an export was done for a study schema that contained the Reporting and Analysis module installed in the same database as the study, messages might appear, saying that the database tables contain a ROWID column, and that values might be obsolete. You can ignore these messages. The messages are related to the materialized views that are being imported. The reporting upgrade script you run later removes the materialized views that were imported.

Note: If an export was done for a study schema that contained the Reporting and Analysis module installed in a different database from the study, an error message might appear, saying that the system could not create a database link. You can ignore this message. This database link is not necessary for the study. If it exists after import, it is removed by the reporting upgrade scripts that you run later.

Rights assigned to right groups

Before upgrading, review and record the rights assigned to each rights group. This information is on the Rights Group detail page in the Admin user interface. For more information, see the *User Guide*. Make sure the rights are correctly assigned for each rights group.

After the upgrade is complete, you use this information to verify that the rights assigned to each rights group are correct. In some cases, the rights assigned to each rights group may not be correct after the upgrade.

Note: This issue can occur only when you upgrade studies from any InForm 4.6 SP2 or earlier release or any InForm 5.0 release. This issue does not occur when you upgrade studies from any InForm 4.6 SP3 or later release, any InForm 5.5 release, or any InForm 6.0 release.

For more information, see *Verify the rights assigned to each rights group* (on page 41).

Itemset types and migrated studies

Two types of itemsets are now available for submitting multiple rows of data on a form at one time: Add Entry itemsets and Repeating Data itemsets.

These itemsets have attributes that are configured during study design. The display and behavior of the itemsets differ in the InForm application depending on the type. Both itemset types can appear on the same form.

Itemsets created in previous InForm releases behave as Add Entry itemsets in the InForm 6.0.0.2 release.

For more information on itemset types, see the *User Guide*.

Imported child controls and migrated studies

In the InForm 6.0.0.2 release, the InForm Data Import utility does not allow the import of dynamic child controls when the parent is not selected.

In releases prior to the InForm 6.0 release, you could use the InForm Data Import utility to import child controls when the parent control was not selected. As a result, if you migrate a study with imported child control data into the InForm 6.0.0.2 release, when a form loads, the dynamic control functionality correctly hides child controls with unselected parent controls, but it does not clear the data in the child controls.

For more information on dynamic control functionality, see the *User Guide*.

Changed resource files

If you customized any resource files, back them up before you upgrade to avoid losing your customizations.

After upgrading, you should review any customizations you made to resources in your previous release and determine whether to reapply them to this release. To reinstall customized resource files, use the MedML Installer utility.

Files updated in InForm release 6.0.0.2

...\XMLBase

sysVolume.xml

Help/Version.htm

Files updated in InForm release 6.0

...\Resource_UI

AdvadedFilterSupport.html

Calendar.js

ContextMenuFrame.html

ControlPanelNoFrame.html

CRFFrame.html

DetailFrame.html

EQACFrame.html

EQACQueryFrame.html

FORGOTPWD.HTML

FORGOTPWDCORRECT.HTML

formhelper.js

Login.css

Login.html

mainstyle.css

MonitorFrame.html

MultiPage.html

Pfscript.js

QueryListFrame.html

Reauth.html

Resetpwdnotification.txt
Resource_auth_template.xml
Resources_template.xml
ScreeningListFrame.html
ScreeningLog.html
SDVListTileBody.html
SDVTitleNoFrame.html
SignPage.html
SSORedirect.html
TabControlframe.html
Top_Reauth.js
TopFrameset.html
TopWindowClose.html
VisitBar.css
VisitBarScript.js
Visitcalculator.html
...\XMLBase
Clintrial_Users.xml
core_ReprtingUser.xml
core_Rights.xml
core_RightsGroups.xml
core_SystemConfig.xml
core_SystemConfigThing.xml
core_Users.xml
core_UsersType.xml
Default_RightsGroups.xml
Defaults.rsp
syscore.rsp
sysform_AdminUserNew.xml
sysform_DefaultConnectionSettings.xml
sysform_Password.xml
sysform_RightsEdit.xml

sysform_SystemConfig.xml
sysform_SystemInfo.xml
sysform_UserType.xml
sysVolume.xml
sysVolume.xmltmpl
upgradeAdminConfigAttributes.xml
UpgradeTrialResources.rsp

New files in InForm release 6.0

...\Resource_UI

ExternalLoginWindowClose.html
ItemsetRowCommentFormFrameset.html
ItemsetRowCommentInterFrameset.html
ReauthenticationWindow.html
Top_Reauth_Trust.js

...\XMLBase

sysform_ItemsetRowComment.xml
UAUM_Users.xml
UpgradePFReportingUser.xml

Files removed in InForm release 6.0

...\Resource_UI

ArchitectCRFFrameset.html
ArchitectCRFFrameset2.html
ArchitectCRFTitle.html
ArchitectCRFTitleRepeating.html
ArchitectFrameset.html
ArchitectItemSetFrameset.html
GenericFrame2Panels.html
ListingsTitle.html
LIVECHECKBOX.HTMl
LIVERADIOBUTTON.HTMl
NOEDITCHECKBOX.HTMl

NOEDITRADIOBUTTON.HTML
PRINTCHECKBOX.HTML
PRINTRADIOBUTTON.HTML
VERTDEEPLIVECHECKBOX.HTML
VERTDEEPLIVERADIOBUTTON.HTML
VERTDEEPNOEDITCHECKBOX.HTML
VERTDEEPNOEDITRADIOBUTTON.HTML
VERTDEEPPRINTCHECKBOX.HTML
VERTDEEPPRINTRADIOBUTTON.HTML
VERTLIVECHECKBOX.HTML
VERTLIVERADIOBUTTON.HTML
VERTNOEDITCHECKBOX.HTML
VERTNOEDITRADIOBUTTON.HTML
VERTPRINTCHECKBOX.HTML
VERTPRINTRADIOBUTTON.HTML

...\XMLBase

core_MercRightsGroups.xml
core_MercUsers.xml
Default_Users.xml

Create backup files

- 1 Stop all studies and servers, using the following command:

PFADMIN STOP SERVER <servername> /Trials

- 2 Back up all customized screens, images, reports, and database scripts, such as views, so they can be reapplied after the upgrade.

Also, back up any other files you might have created or copied to the InForm folder tree, including files in the **InForm\Trials** folder.

Examples:

- Changes to reason drop-down lists on CRF detail pages.
- Logo files and other images.

For more information, see *Database and file customizations* (on page 11).

- 3 If you are using the Reporting and Analysis module or if you made any customizations, you must record or export the registry values.
 - For InForm 4.6, 5.0, or 5.5 releases, record the current the registry values in the **HKEY_LOCAL_MACHINE\SOFTWARE\Phase Forward\AuthenticationFilter** key so you can edit the registry on the new server to apply your customized settings.
 - For InForm 6.0 releases, export the registry values in the **HKEY_LOCAL_MACHINE\SOFTWARE\OracleHS\AuthenticationFilter** key so you can import the values into the registry on the new server.

You will use this information when you restore the Windows registry settings, which is part of the migration process. For more information, see *Restore the Windows registry settings* (on page 37).

- 4 Back up all of the files for each study.

These objects are located in the <InstallationDirectory>\InForm\Trials folder.

- 5 If you are upgrading from an InForm 4.6, 5.0, or 5.5 release, back up the study CSML file located in the <InstallationDirectory>\InForm\Trials folder.

You can use this file to create a new deployment package in the Central Designer 2.0 application. The CSML file uses this naming convention:

<StudyName>_<StudyVersion>.csml (for example, pfst55_0.0.2.csml).

For more information, see *Create a Central Designer 2.0 deployment package* (on page 23).

- 6 If you are upgrading from an InForm 6.0 release, create a backup of the Central Designer rule engines.

The rule engine for studies that are created in the Central Designer application is saved in the <InstallationDirectory>\InForm\Bin\DesignerVersion<Version_Number> folder.

- To create a backup, copy all the files in every DesignerVersion<Version_Number> folder to another location.

These folders include the Central Designer rule engine, external functions, COM interoperability libraries, and configuration files.

- 7 If the CIS software is installed on the InForm application server, you must uninstall it.

Use the **Control Panel > Add/Remove Programs** option to remove the CIS software.

Prepare the InForm 4.6 study database schema for backup

If you are migrating an InForm 4.6 study, you must prepare the InForm 4.6 study database schema for backup before migrating to the InForm 6.0.0.2 release.

Note: Character set differences between the InForm 4.6 release and the InForm 6.0.0.2 release might cause issues with table column widths. Perform this step to prevent these issues.

- 1 Copy the InForm\InstallSupport folder from the downloaded ISO image to the InForm 4.6 application server.
- 2 Run **UpgradeTrialSchema.cmd** (located in the InForm\InstallSupport folder) on the InForm 4.6 application server to prepare the InForm 4.6 study schemas for backup, using the following syntax:

```
UpgradeTrialSchema.cmd <TNS_Service_Name>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **orasystem_user**—Oracle database system user.
- **orasystem_user_pass**—Password for the Oracle database system user.
- **trial_user**—Owner of the study database schema.
- **trial_user_pass**—Password for the owner of the study database schema.

You can also pass the parameters by using the **/accountparams:***"path_to_parameter_file"* command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

The UpgradeTrialSchema command creates a log file **Install_upgrade_log_<current-date>_<current-time>.txt** in the current directory.

Note: If error messages indicating that the system cannot find the paths specified for the System Product Locale and the System Study Locale appear, disregard the messages and continue with the migration.

Back up the InForm study database

Back up the InForm 4.6 or 5.0 database and customized files—Oracle 10g

If you are migrating an InForm 4.6 or InForm 5.0 study from an Oracle 10g database, backing up the InForm study database is required. You will import the DMP file you create in this step during the migration to the InForm 6.0.0.2 release. For more information, see *Import the DMP files* (on page 27).

To back up an InForm 4.6 or InForm 5.0 study database to a DMP file:

- 1 Copy the InForm\InstallSupport folder from the downloaded ISO image to the InForm application server.
- 2 Run **ExportDBUpgrade.cmd** (located in the InForm\InstallSupport folder) on the InForm application server to back up each InForm study database as a DMP file, using the following syntax:

```
ExportDBUpgrade.cmd <DMP_file> <TNS_Service_Name>
```

where:

- **DMP_file**—Name of the DMP file, using the full path.

Note: Do not include spaces in the DMP file name.

- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **trial_user**—Owner of the InForm study database schema.
- **trial_user_password**—Password for the owner of the InForm study database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

Example:

```
ExportDBUpgrade.cmd D:\PFST46.dmp TRIAL1
```

The ExportDBUpgrade command creates a log file **<DMP_file>.log** in the current directory.

Back up the InForm 4.6, 5.5, or 6.0 database and customized files—Oracle 11g

If you are migrating an InForm 4.6, 5.5, or 6.0 study from an Oracle 11g database, backing up the InForm study database is required. You will import the DMP file you create in this step during the migration to the InForm 6.0.0.2 release. For more information, see *Import the DMP files* (on page 27).

To back up an InForm 4.6, 5.5, or 6.0 study database to a DMP file:

- 1 Copy the InForm\InstallSupport folder from the downloaded ISO image to the InForm database server.
- 2 Run **ExportDB.cmd** (located in the InForm\InstallSupport folder) on the InForm database server to back up each InForm study database as a DMP file, using the following syntax:

```
ExportDB.cmd <TNS_Service_Name> <Physical_Directory_Path> <DMP_File_Name>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.
- **Physical_Directory_Path**—Location of the directory that you created on the InForm database server for the data pump.

Note: Do not include spaces in the directory path to the DMP file.

- **DMP_File_Name**—Name of the DMP file to export.

Note: Do not include spaces in the DMP file name.

You are prompted for the following values:

- **orasystem_user**—Name of the Oracle database system user.
- **orasystem_user_pass**—Password for the Oracle database system user.
- **trial_user**—Name of the study database schema owner.
- **trial_user_pass**—Password for the owner of the study database schema.

You can also pass the parameters by using the `/accountparams:"path_to_parameter_file"` command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is `parameter=value`. There is a new line for each parameter, and there are no spaces on a line.

For example:

```
ExportDB.cmd trial11 /u01/app/oracle/dp_export pfst55.dmp
```

The ExportDB command creates the log file `<DMP_File>.log`.

Record study, customer-defined database, and connection information

- 1 Run the pfadmin command using the following syntax:
`pfadmin view service`
- 2 Record the following information for each server and study:
 - Study name, UID and PID.
 - CDD DSN, UID and PID.
- 3 Record the database connection information.

Create a Central Designer 2.0 deployment package

Note: This step is not required if you are upgrading from any InForm 6.0 release.

The InForm 6.0.0.2 software is a 64-bit application. It cannot use the existing 32-bit Central Designer rules engine from any InForm 4.6 5.0, or 5.5 release. If you are upgrading from any of these releases, you must create a new deployment package for each study in the Central Designer 2.0 application.

You can create the deployment package two ways:

- From the Central Designer project that corresponds to the deployed version of the study you plan to migrate to the InForm 6.0.0.2 application.
- Create the deployment package from the CSML file for the study you plan to migrate to the InForm 6.0.0.2 application.

The CSML file for a study is located in the

`<Installation_Directory>\InForm\Trials\<StudyName>\<StudyName>_<StudyVersion>.csml`

For more information, see the Central Designer *InForm Design Guide*.

Note: During creation of the package make sure all the languages in the original study are selected.

Migrating InForm studies to a new machine

You can migrate any InForm 4.6, 5.0, 5.5, or 6.0 study to a new machine.

- Oracle does not support an in-place upgrade for any InForm 4.6, 5.0, or 5.5 release.
- For InForm 6.0 releases, you can migrate to a new machine or you can perform an in-place upgrade.

For more information, see *Running an in-place upgrade for any InForm 6.0 release* (on page 38).

Note: Depending on the InForm version you are migrating from, some steps may not be required.

- 1 *Install the InForm 6.0.0.2 software* (on page 25).
- 2 *Create the database user* (on page 26).
- 3 *Create tablespaces other than INFORM* (on page 26).
- 4 Import the DMP files.
 - *Import the InForm 4.6 or InForm 5.0 DMP files—Oracle 10g* (on page 27).
 - *Import the InForm 4.6, 5.5, or 6.0 DMP files—Oracle 11g* (on page 28).
- 5 *Grant the database user privileges* (on page 29).
- 6 *Set up the InForm servers, studies, and customer-defined databases* (on page 29).
- 7 *Upgrade the study schemas to the InForm 6.0.0.2 release* (on page 30).
- 8 *Update the statistics for the InForm studies* (on page 31).
- 9 *Start the InForm server* (on page 31).
- 10 *Create new database objects for the review schema* (on page 31).
- 11 *Set the MEDMLINSTALLER server* (on page 32).
- 12 *Restore the Central Designer rule engines to the original location* (on page 32).

Note: This step is required only if you are upgrading from any InForm 6.0 release.

- 13 *Run the 64-bit Central Designer deployment package* (on page 33).

Note: This step is not required if you are upgrading from any InForm 6.0 release.

- 14 *Upgrade the study resources* (on page 33).

Note: This step is not required if you are upgrading from any InForm 6.0 release.

- 15 *Terminate obsolete users* (on page 34).

Note: This step is not required if you are upgrading from any InForm 6.0 release.

- 16 *Set the locale for InForm 4.6 sites and users* (on page 35).
- 17 *Delete PF_PATH table records* (on page 36).

Note: This step is not required if you are upgrading from any InForm 6.0 release.

- 18 *Perform the final study configuration* (on page 36).
- 19 *Restore the Windows registry settings* (on page 37).
- 20 *Install or upgrade the InForm Adapter software* (on page 37).
- 21 *Start the InForm studies* (on page 37).

Install the InForm 6.0.0.2 software

Note: Before you install the InForm 6.0.0.2 software, make sure your environment meets the requirements for this release. For more information, see the *Release Notes*.

- Install and configure the InForm software, using the Installation Wizard, as described in the *Installation Guide*.

On the Database Configuration page:

- Select the **Install Admin DB** option.
- Select the **Prep Oracle** option.

Create the database user

- Run **DropDB.cmd** (located in the <Installation_Directory>\InForm\Bin\DBOra folder) to create the database user, using following syntax:

DropDB.cmd <TNS_Service_Name>

where:

- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **pfdbauid**—InForm database admin user.
- **pfdbapid**—Password for the InForm database admin user.
- **trial_user**—Owner of the study schema created by the import.
- **trial_user_pass**—Password for the owner of the study schema created by the import.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

For example:

DropDB.cmd TRIAL1

The DropDB command creates the log file **dropdb.err**.

Create tablespaces other than INFORM

Create the tablespaces, other than INFORM, used by the original study, using Oracle tools/commands. The INFORM tablespace exists by default.

For more information on creating the additional tablespaces, see the *Installation Guide*.

Import the DMP files

Import the InForm 4.6 or InForm 5.0 DMP files—Oracle 10g

Before importing the DMP files, see:

- *Database and file customizations* (on page 11).
- *Oracle 11g database import recommendations* (on page 13).

Note: Oracle does not support migrating tables created outside of the standard InForm application schema or by Enterprise Adoption customers. These tables might cause errors when you import the study DMP files.

- If you are migrating an InForm 4.6 or InForm 5.0 study from an Oracle 10g database, run **ImportDBUpgrade.cmd** (located in the <InstallationDirectory>\InForm\Bin\DBOra folder) to import the studies from the study schema DMP files, using the following syntax:

```
ImportDBUpgrade.cmd <DMP_file> <TNS_Service_Name> <From_Trial_UserID>
<To_Trial_UserID>
```

where:

- **DMP_file**—Name of the file containing the study backup.
- **TNS_Service_Name**—TNS name for the database instance.
- **From_Trial_UserID**—Owner of the database schema being migrated.
- **To_Trial_UserID**—Owner of the study schema created by the import.

You are prompted for the following values:

- **orasystem_user**—Oracle database system user.
- **orasystem_user_pass**—Password for the Oracle database system user.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

Example:

```
ImportDBUpgrade.cmd D:\PFST46.dmp TRIAL1 PFST46UID PFST46UID
```

The ImportDBUpgrade command creates a log file **<DMP_file>.log** (for example, pfst46.log).

Note: Disregard import errors regarding the database link and the proxy user. These errors only appear if the study that you have exported included the Reporting and Analysis module. In addition, disregard import errors regarding rule triggers. Rule triggers are refreshed during the study schema update.

Import the InForm 4.6, 5.5, or 6.0 DMP files—Oracle 11g

Before importing the DMP files, see:

- *Database and file customizations* (on page 11).
- *Oracle 11g database import recommendations* (on page 13).

Note: Oracle does not support migrating tables created outside of the standard InForm application schema or by Enterprise Adoption customers. These tables might cause errors when you import the study DMP files.

- If you are migrating an InForm 4.6, 5.5, or 6.0 study from an Oracle 11g database, run **ImportDB.cmd** (located in the `<InstallationDirectory>\InForm\Bin\DBOra` folder) to import the studies from the study schema DMP files, using the following syntax:

```
ImportDB.cmd <TNS_Service_Name> <Physical_Directory_Path> <DMP_File_Name>
<From_User>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.
- **Physical_Directory_Path**—Location of the directory that you created on the InForm database server for the data pump.
- **DMP_File_Name**—Name of the DMP file to import.
- **From_User**—Owner of the database schema being migrated.

You are prompted for the following values:

- **orasystem_user**—Oracle database system user.
- **orasystem_user_pass**—Password for the Oracle database system user.
- **trial_user**—Owner of the study schema created by the import.
- **trial_user_pass**—Password for the owner of the study schema created by the import.

You can also pass the parameters by using the `/accountparams:"path_to_parameter_file"` command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is `parameter=value`. There is a new line for each parameter, and there are no spaces on a line.

For example:

```
ImportDB.cmd trial11 /u01/app/oracle/dp_export pfst46.dmp pfst46uid
```

The ImportDB command creates the log file `<DMP_file>.log`.

Grant the database user privileges

- Run **GrantAccess.sql** (located in the <Installation_Directory>\InForm\Bin\DBOra folder) using the following syntax:

```
sqlplus <pfdbauid>@<trial_instance> @GrantAccess.sql <trial_userid>
```

where:

- **trial_userid**—The name of the database user.

When prompted, enter the InForm database admin user password.

For example:

```
sqlplus infadmin@trial1 @GrantAccess.sql PFST46UID
```

Set up the InForm servers, studies, and customer-defined databases

- 1 To define each server, use the following command and syntax:

```
PFADMIN SETUP SERVER <servername>
```

- 2 To define each study, use the following command and syntax:

```
PFADMIN SETUP TRIAL <trialname> <servername> /db <oracle_connect_string>
```

where:

- **trialname**—Name of the study you are defining.
- **servername**—Name of the server on which the study is located.
- **oracle_connect_string**—String that the InForm server uses to connect to the Oracle Instance for the study.

You are prompted for the following values:

- **uid**—Owner of the study database schema.
- **pid**—Password for the owner of the study database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

- 3 To set up the customer-defined databases using the NoSchema option, use the following command and syntax:

```
PFADMIN SETUP CDD <Refname> <trialname> /db <oracle_connect_string> <DSN>  
[NoSchema]
```

where:

- **RefName**—RefName of the CDD.
- **trialname**—Name of the study you are defining.
- **oracle_connect_string**—String that the InForm server uses to connect to the Oracle Instance for the study.

- **DSN**—Data Source Name of the CDD.
- **[NoSchema]**—Indicates that no new CDD schema should be created during setup.

You are prompted for the following values:

- **uid**—Owner of the study database schema.
- **pid**—Password for the owner of the study database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is **parameter=value**. There is a new line for each parameter, and there are no spaces on a line.

Upgrade the study schemas to the InForm 6.0.0.2 release

The **UpgradeTrialSchema** command upgrades studies from the InForm 4.6 SP2 release and later, all InForm 5.0, 5.5, and 6.0 releases to the InForm 6.0.0.2 release. The command terminates if the InForm release is different from these releases.

- Run **UpgradeTrialSchema.cmd** (located in the *<InstallationDirectory>\InForm\Bin\DBOra* folder) on the InForm 6.0.0.2 application server to upgrade each study schema to the InForm 6.0.0.2 release, using the following syntax:

```
UpgradeTrialSchema.cmd <TNS_Service_Name>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **orasystem_user**—Oracle database system user.
- **orasystem_user_pass**—Password for the Oracle database system user.
- **trial_user**—Owner of the study database schema.
- **trial_user_pass**—Password for the owner of the study database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is **parameter=value**. There is a new line for each parameter, and there are no spaces on a line.

For example:

```
UpgradeTrialSchema.cmd TRIAL1
```

The **UpgradeTrialSchema** command creates a log file **Install_upgrade_log_<current-date>_<current-time>.txt** in the current directory.

Update the statistics for the InForm studies

- Run **updatec.sql** for each study, using the following command and syntax:
`SQLPLUS <pfdb_admin>@<trial_instance> @updatec.sql <trialuserid>`

When prompted, enter the password for the InForm database administrator.

Start the InForm server

- Start the InForm server for the studies, using the following command:

```
PFADMIN START SERVER <servername>
```

Create new database objects for the review schema

Before you create the new database objects for the review schema, see *Database and file customizations* (on page 11).

Note: The review schema requires about two times the space as the study schema. Be sure to check your disk space and tablespace allocations to ensure adequate resources before you create the review schema tables.

Note: The time needed for creating the review schema tables is proportionate to the size of your study. This step can take quite a while.

- Run the **pfadmin RECREATEREVIEWSCHEMA** command to apply any fixes to the review schema and to create new tables for the review schema, using the following syntax:

```
pfadmin RECREATEREVIEWSCHEMA <Trial_Name> <StudyLocale_ISO_name>
```

where:

- Trial_Name**—Name of the study for which you are creating the review schema.
- StudyLocale_ISO_name**—Name of the InForm Study Locale language.

For InForm 4.6 studies, you must use en-US as the *<StudyLocale_ISO_name>*.

For InForm 5.0, 5.5, and 6.0 studies you can pick any of the supported Study Locales.

The pfadmin RECREATEREVIEWSCHEMA command creates the following log files in the *<InstallationDirectory>\InForm\Bin\DBOra\ReviewSchema* folder:

- Install_Clinical_<trialname>_<currentdate>_<currenttime>.txt
- Install_Fixed_<trialname>_<currentdate>_<currenttime>.txt

Note: If any schema objects such as views, functions, procedures, packages, or triggers are invalid, pfadmin RECREATEREVIEWSCHEMA fails. You must repair any invalid schema objects before you can rerun pfadmin RECREATEREVIEWSCHEMA.

Set the MEDMLINSTALLER server

You must set the MEDMLINSTALLER to the current server in order to install the new resources.

The study must be in Design mode before running the commands.

- 1 Start each study in Design mode, using the following command and syntax:

```
PFADMIN START TRIAL <trialname> /DESIGN
```

- 2 Set the MEDMLINSTALLER machine for each study, using the following commands and syntax:

```
PFADMIN SETSERVER MEDMLINSTALLER <trialname> <computername>
```

```
PFADMIN SETSERVER SITE <trialname> <computername>
```

where:

- **computername**—The unqualified name of the host machine.

Restore the Central Designer rule engines to the original location

Note: This step is required only if you are upgrading from any InForm 6.0 release.

- 1 Copy the DesignerVersion<version number> folder or folders you backed up to the InForm installation directory <InstallationDirectory>/InForm/bin.
- 2 Register the Central Designer rule engine DLL for each DesignerVersion<version_number> folder that you copied.
 - a Check on your system for the location of the latest copy of RegAsm.exe, version 2.0.50727, typically located at: **C:\WINDOWS\Microsoft.NET\Framework\v2.0.50727**.
 - b For each folder that you restore, execute the following command, all on one command line:

```
C:\Windows\Microsoft.NET\Framework\v2.0.50727\RegAsm.exe /codebase
<InstallationDirectory>\InForm\bin\DesignerVersion<version
_number>\PhaseForward.Designer.InFormRuleEngine.dll
```

where:

- **InstallationDirectory**—The directory in which the InForm software is installed.
- **DesignerVersion<version_number>**—The name of each DesignerVersion folder that you copied from the bin folder.

Note: The **/codebase** flag is required. If you run RegAsm without the codebase flag, the rule engine registers, but rules do not run properly in the study. This prevents the enrollment of new subjects and the submission of some forms.

- 3 Restore the files contained in the folder <InstallationDirectory>\InForm\Trials\<trial_name>.

Run the 64-bit Central Designer deployment package

Note: This step is not required if you are upgrading from any InForm 6.0 release.

- Run the Central Designer deployment package wizard created in the Central Designer 2.0 application, and supply the following information for the study:
 - InForm server.
 - Study name.
 - Study database schema owner name and password.
 - Study database connection string or study DSN.

Note: You may encounter one or more **Element already exists - elements cannot be changed** errors when you try to deploy an InForm 4.6 study to the InForm 6.0.0.2 application. For a workaround, contact Global Customer Support.

Upgrade the study resources

Note: This step is not required if you are upgrading from any InForm 6.0 release.

- 1 Run **UpgradeResources.cmd** (located in the *<InstallationDirectory>\InForm\Bin\DBOra* folder) to upgrade the resources and the XMLBase changes for each study, using the following syntax:

```
UpgradeResources.cmd <trialname> <InForm_version_upgraded_from> [SILENT]
```

where:

- **InForm_version_upgraded_from**—4.6, 5.0, or 5.5.
 - **[SILENT]**—Option that suppresses parameters that are normally supplied interactively with the command.
- 2 During the upgrade process from the InForm 4.6 release, you are asked if you want to change the System Config QueryMaxLength parameter to 350 characters.
 - Type **Y** and press ENTER to change the value.

If you used the SILENT option when you ran the command, the QueryMaxLength is changed automatically.

The product locale and study locale of the system user are set to en-US if they do not already exist.

Terminate obsolete users

Note: This step is not required if you are upgrading from any InForm 6.0 release.

- Run the following command to terminate built-in system users and remove invalid email addresses:

```
TerminateUsers.cmd <TNS_Service_Name> <trialname>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.
- **trialname**—Name of the study.

You are prompted for the following values:

- **trial_user**—Owner of the study database schema.
- **trial_user_password**—Password for the study database schema.

You can also pass the parameters by using the **/accountparams:***"path_to_parameter_file"* command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is **parameter=value**. There is a new line for each parameter, and there are no spaces on a line.

For example:

```
TerminateUsers.cmd trial1 pfst55
```

Set the locale for InForm 4.6 sites and users

Note: This step is required only when you migrate InForm 4.6 studies to the InForm 6.0.0.2 release.

InForm 5.0, 5.5, and 6.0 studies have locale attributes associated with sites and users, but InForm 4.6 studies do not. When you migrate studies from the InForm 4.6 release to the InForm 6.0.0.2 release, these properties must be added to existing sites and users.

- 1 Run **PFInFormMigrate.cmd** (located in the *<InstallationDirectory>\InForm\Bin* folder) to associate English locale attributes to sites and users, and to create the Upgrade MedML file **UpgradeSitesAndUsers.xml**, using the following syntax:

```
PFInFormMigrate <TNS_Service_Name>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **uid**—Owner of the database schema.
- **pid**—Password for the owner of the database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is **parameter=value**. There is a new line for each parameter, and there are no spaces on a line.

Example:

```
PFInFormMigrate TRIAL1
```

The Upgrade MedML file **UpgradeSitesAndUsers.xml** is created in the *<InstallationDirectory>\InForm\XMLBase* folder.

- 2 Open the **UpgradeSitesAndUsers.xml** file you created in step 1, and check the site MedML definitions for warnings on the **TIMEZONE** attribute.

If the Time Zone specified for a site in the InForm 4.6 study does not match the time zones defined for the operating system on the InForm application server, the MedML for the given site contains a warning for the **TIMEZONE** attribute along with the original Time Zone.

- 3 If a warning appears, change the **TIMEZONE** attribute to a valid Time Zone key name in the InForm server registry before proceeding to the next step.

Time Zones are located at

HKEY_LOCAL_MACHINE\SOFTWARE\Microsoft\Windows NT\CurrentVersion\Time Zones in the registry.

- 4 Ensure the XML file is well formed after the changes.

- Run **UpgradeSitesAndUsers.cmd** (located in the <InstallationDirectory>\InForm\Bin\DBOra folder) to install UpgradeSitesAndUsers.xml, using the following syntax:

```
UpgradeSitesAndUsers.cmd <TrialName>
```

where:

- **TrialName**—The name of the study being migrated.

For example, run the following command to upgrade sample study PFST46:

```
UpgradeSitesAndUsers.cmd PFST46
```

The UpgradeSitesAndUsers command sets the product locale and study locale of all users except the system user to en-US. Similarly, all sites have site locale set to en-US and the User Name Order property set to **F** (given name), **L** (surname).

Delete PF_PATH table records

Note: This step is not required if you are upgrading from any InForm 6.0 release.

For an upgraded study, the PF_PATH table in the study schema may include records that are not required by the InForm 6.0.0.2 software. These records may cause unpredictable results. The following commands delete all the records from the PF_Path table so they can be repopulated when the study restarts:

- Stop the study.

```
pfadmin stop trial <trialname>
```

- Run the clearcontrolpaths.sql script (located in the <Installation_Directory>\InForm\Bun\DBOra folder) using the following syntax:

```
sqlplus <trial_userid>@<TNS_Service_Name>@clearcontrolpaths.sql
```

When prompted, enter the password for the study user.

where:

- **trial_userid**—Owner of the study database schema.
- **TNS_Service_Name**—TNS name for the database instance.

Example:

```
sqlplus pfst60uid@trial1@clearcontrolpaths.sql
```

Perform the final study configuration

- Finish configuration of the studies related to customer-defined database and randomization.
For more information, see the *Study and Reporting Setup Guide*.
- Reapply any customizations.
For more information, see the *Study and Reporting Setup Guide*.

Restore the Windows registry settings

- 1 Restore the Windows registry settings that you recorded or exported when you created the backup files to the new registry key
HKEY_LOCAL_MACHINE\SOFTWARE\OracleHS\AuthenticationFilter.

For more information, see *Create backup files* (on page 19).

- For InForm 4.6, 5.0. or 5.5 releases, add the registry values that you recorded from the **HKEY_LOCAL_MACHINE\SOFTWARE\Phase Forward\AuthenticationFilter** key.
- For InForm 6.0 releases, import the registry values that you exported from the **HKEY_LOCAL_MACHINE\SOFTWARE\OracleHS\AuthenticationFilter** key.

- 2 Restart IIS.

Install or upgrade the InForm Adapter software

Perform this step if you are currently using, or plan to start using, the InForm application with products that require the InForm Adapter software.

To install or upgrade the InForm Adapter software, refer to the InForm Adapter *Installation and Configuration Guide* and the InForm Adapter *Release Notes*.

For information about release compatibility, see the *Release Notes*.

Note: The InForm Adapter software and the InForm Server Adapter software, as well as the documentation images, are available from the Download Center.

Start the InForm studies

- Start each study, using the following command and syntax:

```
PFADMIN START TRIAL <trialname>
```

Running an in-place upgrade of any InForm 6.0 release

Use the following steps to run an in-place upgrade of any InForm 6.0 studies and the reporting environment.

- 1 Download the InForm 6.0.0.2 ISO file and extract the contents to a location on the InForm application server.
- 2 Open a Command Prompt window.
- 3 Stop the studies and servers, by typing the following commands in the order listed.


```
PFADMIN STOP TRIAL <studyname>
```

```
PFADMIN STOP SERVER <servername>
```
- 4 Stop all InForm applications and integrations that require the InForm service, for example, the Central Coding application.
- 5 Stop the InForm service.
- 6 Back up the following:
 - The contents of the <Installation_Directory>\InForm folder.
 - All of the files for each study.

These objects are located in the <Installation_Directory>\InForm\Trials folder.
- 7 Export the current registry settings from all servers, including the InForm server and the Reporting and Analysis server in **HKEY_LOCAL_MACHINE\SOFTWARE\OracleHS**.

Note: Some of these values are reset to their defaults when you re-install the InForm software or when you run the in-place upgrade.

- 8 Navigate to the folder where you extracted the contents of the ISO file, and) run the **setup.exe** file located in the InFormUpgrade folder.

The installation wizard opens.
- 9 Follow the prompts to run the wizard.

Running the wizard upgrades the core InForm software.

 - If you are not using Reporting, the upgrade is complete. Proceed to Step 11.
 - If you are using Reporting, proceed to Step 10.
- 10 Copy **Op package and reports 6_0_0_2.zip**, located in the <Installation_Directory>\InForm\Bin\DBOra\Reporting folder on the InForm application server, to the <Installation_Directory>\deployment folder on the Cognos server.

Note: This step upgrades your Reporting environment.

- 11 Restore the HKEY_LOCAL_MACHINE\SOFTWARE\OracleHS registry settings from the file you exported in Step 7.
- 12 Start the InForm service.
- 13 Start the installed studies and servers, by typing the following commands in the order listed.

```
PFADMIN START SERVER <servername>
```

```
PFADMIN START TRIAL <studyname>
```

Post-upgrade tasks

Clear the client-side browser caches

Clear the Internet Explorer 7.0 browser cache

- 1 Select **Tools > Internet Options**.
- 2 In the **Browsing history** section of the **General** tab, click **Delete**.
The Delete Browsing History dialog box appears.
- 3 Click **Delete Files**.
A confirmation dialog box appears.
- 4 Click **Yes**.
- 5 Click **Delete Cookies**.
A confirmation dialog box appears.
- 6 Click **Yes**.
- 7 Click **Close**.
- 8 Click **OK**.

Clear the Internet Explorer 8.0 or 9.0 browser cache

- 1 Select **Tools > Internet Options**.
- 2 In the **Browsing history** section of the **General** tab, click **Delete**.
- 3 Deselect **Preserve Favorites website data**.
- 4 Select **Temporary Internet Files**.
- 5 Select **Cookies**.
- 6 Click **Delete**.
- 7 Click **OK**.

Clear the Apple Safari 5.1 browser cache

- 1 Select **Safari > Empty cache**.
- 2 Click **Empty**.

Clear the Firefox 3.6 or 11.0 browser cache

- 1 Select **Tools > Clear Recent History**.
The Time range to clear window appears.
- 2 From the drop-down list, select **Everything**.
- 3 In the **Details** list, select **Cache**.
- 4 Click **Clear Now**.

Clear the Chrome 19 browser cache

- 1 Select **Settings >Tools > Clear browsing data**.
Clear browsing data window appears
- 2 From the Clear the following items from: drop-down list, select **clear from beginning of time**.
- 3 Select all the items in the list.
- 4 Click **Clear browsing data**.

Verify the rights assigned to each rights group

In some cases, the rights assigned to each rights group may not be correct after the upgrade.

Note: This issue can occur only when you upgrade studies from any InForm 4.6 SP2 or earlier release or any InForm 5.0 release. This issue does not occur when you upgrade studies from any InForm 4.6 SP3 or later release, any InForm 5.5 release, or any InForm 6.0 release.

- 1 After the upgrade is complete, verify that the rights assigned to each rights group are correct using the information you recorded before you upgraded your study.
This information is on the Rights Group detail page in the Admin user interface. For more information, see ***Rights assigned to right groups*** (on page 13) and the *User Guide*.
- 2 Update the rights for each rights group, as needed.

Reset your password

When you log in to the InForm application for the first time after upgrading to the InForm 6.0.0.2 release you are required to reset your password so that you can select the security questions required by the password reset functionality.

Note: This step is not required if you are upgrading from the InForm 4.6.4 release or any InForm 6.0 release.

To log in to the InForm application:

- 1 Open an Internet browser.
- 2 Type the URL for your study, and press Enter.

The Change Expired Password page appears. In addition to the user name, password, and email address fields, this page displays:

- Two drop-down lists from which you select a security question.
There are four pre-defined security questions to choose from. You cannot modify or change the questions.
 - What is your mother's maiden name?
 - What is the name of your pet?
 - What is the city of your birth?
 - What is your favorite color?
- Two text boxes in which you enter the answers for the selected security questions.
You must provide answers to both questions. The answers are case-insensitive. Leading and trailing spaces are ignored.

For more information, see the *User Guide*.

CHAPTER 4

Upgrading and migrating the reporting environment to this release

In this chapter

Overview of upgrading and migrating the reporting environment to this release	44
Migrating the reporting environment to a new machine.....	45
Running an in-place upgrade of an InForm 6.0 reporting environment	59

Overview of upgrading and migrating the reporting environment to this release

You can upgrade the reporting environment for any InForm 4.6 SP2 studies or later, any InForm 5.0, 5.5, or 6.0 studies. For more information, see ***Requirements for the InForm studies you can migrate*** (on page 3).

Upgrading the reporting environment is only required if you are using the Reporting and Analysis module.

For more information, see:

- ***Migrating the reporting environment to a new machine*** (on page 45).
- ***Running an in-place upgrade of an InForm 6.0 reporting environment*** (on page 59).

Migrating the reporting environment to a new machine

You can migrate the reporting environment for any InForm 4.6, 5.0, 5.5, or 6.0 study to a new machine.

- Oracle does not support an in-place upgrade for any InForm 4.6, InForm 5.0, or InForm 5.5 release.

You must install a new Cognos 8.4.1 64-bit environment, and create, export, and then import your export packages.

- For InForm 6.0 releases, you can migrate to a new machine or you can perform an in-place upgrade.

For more information, see *Running an in-place upgrade of an InForm 6.0 reporting environment* (on page 59).

Note: Depending on the InForm version you are migrating from, some steps may not be required.

- 1 *Install and configure the Cognos 8.4.1 environment* (on page 46).
 - a *Install and customize the Cognos 8.4.1 software* (on page 46).
 - b *Identify or create the required users for any InForm 4.6, InForm 5.0, or InForm 5.5 study* (on page 47).
- 2 *Create an export package and export the data from your Cognos environment* (on page 48).
 - a Review *About the export package* (on page 48).
 - b *Create an export package for each study in the Cognos 8.2 environment or the Cognos 8.4 environment* (on page 48).
- 3 *Create the streams between the InForm study schema and the Reporting schema, and install the Reporting schema* (on page 50).
- 4 *Import the export package* (on page 51).
- 5 *Run the PFCognosConfig utility* (on page 55).
- 6 *Modify the PFRSetupTrial.xml file* (on page 55).
- 7 *Run PFRInit on the InForm application server* (on page 56).
- 8 *Delete obsolete reports* (on page 56).

Note: This step is required only if you are upgrading from any InForm 4.6 studies.

- 9 *Associate custom reports with the new operational packages* (on page 57).

Note: This step is required only if you are upgrading from any InForm 4.6 studies.

- 10 *Restart the InForm Model Updater service* (on page 57).
- 11 *Revalidate custom reports* (on page 58).

- 12 *Update schedules and jobs* (on page 58).
- 13 *Delete the study organizational units from the LDAP configuration* (on page 58).

Note: This step is not required if you are upgrading from any InForm 6.0 release.

Install and configure the Cognos 8.4.1 environment

Install and customize the Cognos 8.4.1 software

- 1 Install the Cognos 8.4.1 64-bit software on the Cognos servers.

Note: For information on how to install and customize the Cognos 8.4.1 software, see the *Installation Guide*.

- 2 Install Fix Pack 4 (FP4) for the Cognos 8.4.1 software.
- 3 Install the Cognos hotspots in this order:
 - a up_c8bisrvr_win64_8.4.105.1076
 - b up_c8bisrvr_win64_8.4.105.1100
 - c up_c8bisrvr_win64_8.4.107.1013
- 4 Customize the Cognos 8.4.1 software on the Cognos servers.
 - a Run the Cognos 8 Customization for InForm wizard (CRNConfig) on all Report and Content Manager servers.
 - b Run the Cognos 8 Gateway Customization for InForm wizard (CRNGatewayConfig) on all Gateway servers.
 - c Start the Cognos 8 Business Intelligence servers.

Identify or create the required users for any InForm 4.6, 5.0, 5.5, or 6.0 study

Before you start the Cognos 8.4.1 upgrade, make sure that the following users exist so that you can log in to the Cognos 8 Business Intelligence server:

- A Cognos user who is a System Administrator.
- An InForm study user who is assigned to a rights group that includes the Reports right and is also a member of the Directory Administrators reporting group.
- An InForm study user (pfreportinguser) who is a member of the Sponsor Users and Publishers reporting group.

The pfreportinguser will be used for the first time by the pfmodelupdater service. The pfmodelupdater service uses the pfreportinguser user name and password to generate the clinical model.

- If necessary, configure the pfreportinguser account.
For more information, see the *Installation Guide*.
- Change the password for the pfreportinguser using both the user interface and the pfadmin command.
 - Log in to the InForm application and change the password for the pfreportinguser using the user interface.
 - Change the password for the pfreportinguser using the following command:

```
PFADMIN SETSERVER PFREPORTINGUSERPW <studyname>
```

When prompted, enter the password for the pfreportinguser.

Note: Make sure that the password you enter through the user interface and the command line match.

- Reset the password for the Reporting CAP Admin database using the following command:

```
PFADMIN SETSERVER REPORTINGCAPADMINPW <studyname>
```

When prompted, enter the password for the reportingcapadmin user.

Create an export package and export the data from your Cognos environment

About the export package

You must create reporting export packages for each InForm study being migrated to the InForm 6.0.0.2 release while your current InForm and Cognos environments are still active. The export packages contain:

- A study folder.
- An operational package (InForm Trial Management for <trialname>).
- A clinical package (<trialname> Clinical) from the Cognos Public Folders.

You will import these packages into the Cognos 8.4 environment after you install the InForm 6.0.0.2 software and the Cognos 8.4 software.

Because the current Cognos environment remains intact until the Cognos 8.4 environment is active, there is no need to back up any files on the Cognos server except for the export packages.

Create an export package for each study in the Cognos 8.2 environment or the Cognos 8.4 environment

For each InForm study, ensure that all users move the content in their **My Folders** area to the study folder in **Public Folders**.

Note: Content in the **My Folders** area and individual User Preferences are not migrated during the export and import steps of the Cognos upgrade.

Note: In multi-server Cognos environments, the deployment archives are created on the active Content Manager server.

- 1 Log in to an InForm study as a user who has the Reports right and is a member of the Directory Administrators reporting group.
- 2 Click **Reports**.
- 3 Select **Tools > Logon** to log in to the Cognos 8 Business Intelligence 8.2 or 8.4 server as a System Administrator.
- 4 In the **User** field, type the user name for the Cognos System Administrator user (for example, crnsysadmin).
- 5 In the **Password** field, type the password for the Cognos System Administrator user, and click **OK**.
- 6 For Cognos 8.2:
 - Select **Tools > Content Administration**.For Cognos 8.4:
 - a Select **Launch > Reporting Administration**.
 - b On the Configuration tab > **Content Administration**.

- 7 Click **New Export**.
The Export wizard appears.
- 8 Enter a name for the export package, and complete the fields on the **Description** page.
- 9 Click **Next**.
- 10 On the **Choose a deployment method** page, ensure that the **Select public folder and directory content** option is selected, and click **Next**.
- 11 On the **Select the public folders content** page, in the Public Folders content section, select the **trial folder**, the **operational package** (InForm Trial Management for <trial_name>), and the **clinical package** (<trial_name> Clinical).

Note: If you are using Cognos 8.4, Select Add to before you can select the packages you want to include.

- 12 Update the settings on each page for the export deployment package as listed in the following table.

Section	Option	Setting
Public folders content page		
Public folders content section	Disable after import	Selected for the study folder, the operational package, and the clinical package
Options section	Include report output version	Selected
	Include run history	Selected
	Include schedules	Selected
Directory content page		
Directory content	Include Cognos groups and roles	Deselected
	Include distribution lists and connections	Deselected
	Include data sources and connections	Deselected
General options page		
Access permissions	Include access permissions	Deselected
External namespaces	Do not include references to external namespaces	Selected
Entry ownership	Set the owner to	The user performing the import
	Apply to	New entries only

- 13 Click **Next**.
- 14 On the **Specify a deployment archive** page, click **Next**.
- 15 On the **Review the Summary** page, review the export details, and click **Next**.
- 16 On the **Select an action** page, select **Save and Run once**, and click **Finish**.
- 17 On the **Run with options** page, in the **Time** section, select **Now**, and click **Run**.
- 18 On the next page, after closing the dialog box, select **View the details of this export**, and click **OK**.
The View an export deployment record page appears, and displays the status of the export.
- 19 Click **Close**.

Create the streams between the InForm study schema and the Reporting schema, and install the Reporting schema

Note: You only need to perform this step if your InForm study and the Reporting and Analysis module are installed in different database instances.

- 1 Configure **configdiffdb.sql**.
For more information, see the *Installation Guide*.
- 2 Re-apply any customizations.
 - a If you customized any files using the **streams_table_adds** variable, or by modifying **tabexcludelist.sql** or **tabincludelist.sql**, set the **streams_table_adds** variable in **configdiffdb.sql**.
If these tables already have a prefix of PF*, CT*, or RT*, they will automatically be replicated.
 - b To exclude tables with a prefix of PF*, CT*, or RT* from replication, set the **streams_table_exclude** variable in **configdiffdb.sql**.
- 3 Open a Command Prompt window.
- 4 Navigate to the folder in which the Reporting and Analysis software is located (for example, e:\InForm\Bin\DBOra\Reporting).
- 5 Run **configandcheckdiffdb**.
For more information, see the *Installation Guide*.
- 6 Run **install_reporting_diffdb** to install the reporting schema.
For more information, see the *Installation Guide*.

Import the export package

The Reporting Import utility is a command-line script that imports report content into a new or upgraded environment, or imports custom report specifications into an existing environment. Using this utility eliminates a manual step in the upgrade and migration process, and provides additional security and accuracy.

Use this command to migrate your export package into the InForm 6.0.0.2 environment.

Running the Reporting Import Utility

Run the ImportUtility command (located in the *<Installation_Directory>\InForm\bin* on the InForm application server) using the following syntax:

```
ImportUtility -all
```

When prompted, enter the:

- Cognos system administrator user namespace.
- User name for the Cognos system administrator.
- Cognos system administrator password.
- Internal URI that the InForm server uses to communicate with the Cognos 8 BI server.

This Cognos 8 BI parameter is set in the InForm CRN wizard. The parameter setting is stored in the cogstartup.xml file and corresponds to the Reporting internal URI value on the Admin > System Configuration page of the InForm application.

- Name of the deployment package containing the reporting content you want to import into the Reporting server.
- Password for the deployment package.

Note: You can also pass the system administrator namespace, user name, dispatcher URL, and deployment package name in a parameter file. For more information, see *ImportUtility* (on page 52).

The ImportUtility command:

- Creates a secure staging folder, accessible to users with system administration privileges.
The name of the folder is SecureImport_<timestamp> where the time stamp format is yyyyymmdd_hhmmss_sss, where sss represents milliseconds (for example SecureImport_20130920_172055_121).
- Creates a new import specification object, accessible to users with system administration privileges.
- Imports the content of the deployment package to the SecureImport_<timestamp> folder, and makes it accessible only to users with system administration privileges.
- Copies the imported reporting content to the Reporting and Analysis Public folder.
- Deletes the SecureImport_<timestamp> folder, any content it contains, and the import specification object.

For more information, see *ImportUtility* (on page 52).

ImportUtility

Purpose

Securely import reporting content you export using Cognos commands to a new instance of the InForm Reporting server.

Location

The InForm application server folder on the InForm application server. For example, <Installation_Directory>\InForm\bin.

Usage

ImportUtility [-all | -import | -nopwd | -pfmtr] ["path_to_parameter_file"]

Option	Parameter
-all	<ul style="list-style-type: none">• Creates a secure staging folder on the Reporting server accessible to system administrators. The name of the staging folder is SecureImport_<timestamp>.• Creates the import specification object on the Reporting server accessible to system administrators.• Imports the content from a password-protected deployment package to the SecureImport_<timestamp> folder and makes it accessible only to system administrators.• Copies the imported content to the Public folders on the Reporting server.• Deletes the SecureImport_<timestamp> folder, all its content, and the import specification object.

Option	Parameter
-import	<ul style="list-style-type: none"> Creates a secure staging folder on the Reporting server accessible to system administrators. The name of the staging folder is SecureImport_<timestamp>. Creates the import specification object on the Reporting server accessible to system administrators. Imports the content from a password-protected deployment package to the SecureImport_<timestamp>. Deletes the import specification object.
-nopwd	<ul style="list-style-type: none"> Creates a secure staging folder on the Reporting server accessible to system administrators. The name of the staging folder is SecureImport_<timestamp>. Creates the import specification object on the Reporting server accessible to system administrators. Imports the deployment package content to the SecureImport_<timestamp> folder and makes it accessible only to system administrators. Copies the imported content to the Public folders on the Reporting server. Deletes the SecureImport_<timestamp> folder, all its content, and the import specification object.
-pfmtr	<ul style="list-style-type: none"> Creates a secure staging folder on the Reporting server accessible to system administrators. The name of the staging folder is SecureImport_<timestamp>. Creates the import specification object on the Reporting server accessible to system administrators. Imports the content to the SecureImport_<timestamp> folder. Deletes the import specification object.
path_to_parameter_file	<ul style="list-style-type: none"> When specified, includes the path to a text file that contains the parameters required to run the command. You will always be prompted to enter required passwords on the command line. If the parameter file is not specified, the command prompts for the required parameters. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

Command line prompts

Depending on the option you specify with the ImportUtility, you are prompted for following parameters:

- **sysadmin_namespace**—Cognos system administrator user namespace.
- **sysadmin_uid**—User name for the Cognos system administrator.
- **sysadmin_pass**—Cognos system administrator password.
- **dispatcher_url**—Internal URI that the InForm server uses to communicate with the Cognos 8 BI server.

This Cognos 8 BI parameter is set in the InForm CRN wizard. The parameter setting is stored in the cogstartup.xml file and corresponds to the Reporting internal URI value on the Admin > System Configuration page of the InForm application. For example, <http://example.com:9300/p2pd/servlet/dispatch>.

- **archive_name**—Name of the deployment package containing the reporting content you want to import into the Reporting server.
- **archive_pass**—Password for the deployment package.

Note: The **archive_pass** parameter is not required for the **-nopwd** or the **-pfmtr** option.

Parameter file contents

The following parameters may be passed to the ImportUtility command in a parameter file.

- **sysadmin_namespace**=Cognos system administrator user namespace.
- **sysadmin_uid**=User name for the Cognos system administrator.
- **dispatcher_url**=Internal URI that the InForm server uses to communicate with the Cognos 8 BI server.
- **archive_name**=Name of the deployment package containing the reporting content you want to import into the Reporting server.

Notes

- The **-import** and **-pfmtr** command options should be used to stage custom reports before distributing them to multiple studies.

To move the content from the staging folder to the study folders and set study-specific permissions, use the **pfmtrsetuputil** command.

- The **-nopwd** and **-pfmtr** command options are used when you are importing a deployment package that is not password protected.

These options should be used only when there is no clinical data present in the deployment package.

Example

```
ImportUtility -all "E:\scripts\import.txt"
```

Run the PFCognosConfig utility

Run the PFCognosConfig utility (located in the `<InstallationDirectory>\InForm\bin` folder) on the InForm 6.0.0.2 application server to configure an InForm study to work with the Cognos 8 Business Intelligence application. For more information, see the *Installation Guide*.

Modify the PFRSetupTrial.xml file

The PFRSetupTrial.xml file is located in the `<InstallationDirectory>\InForm\Bin` folder.

- 1 In the **TrialFolder** section, change the **UserRoot** tag to **ReportsUpgrade**.
ReportsUpgrade is a tag that triggers the reporting specification upgrade.
- 2 Change the value of the **connectionstring** from *dbconnectionstring* to the Oracle connection string for the Reporting and Analysis database.
In a single database configuration, use the study database connection string.
- 3 Update the DBUID parameter with the UserID of the study owner.
- 4 When you run the PFRInit command you have to pass the following parameters on the command line:
 - **sysadmin_namespace**—Oracle Directory Server admin namespace.
 - **sysadmin_uid**—User name for the Oracle Directory Server admin namespace.
 - **sysadmin_pass**—Password you chose when creating the Oracle Directory Server admin namespace user.
 - **publisher_namespace**—Custom Authentication Provider (CAP) namespace.
 - **publisher_uid**—InForm study user who is a member of the Publisher and Sponsor Users (or Site Users) Reporting groups.
 - **publisher_pass**—Password of the InForm study user.
 - **dispatcher_url**—Internal URI that the InForm server uses to communicate with the Cognos 8 BI server.
 - **trial_name**—Name of the InForm study.
 - **reportdb_pass**—Password for the reporting schema user.

There is also an option to pass these values in a parameter file. For more information, see the *Installation Guide*.

Note: You must use the most current version of this file, located in the `<InstallationDirectory>\InForm\bin` directory for this release. Do not use a backed up version of this file from an earlier release of the InForm software.

Run PFRInit on the InForm application server

- Run the PFRInit command from the <InstallationDirectory>\InForm\Bin folder.

For more information, see the *Installation Guide*.

The PFRInit command:

- Creates a new operational package called <studyname> **InForm Trial Management** in the Public Folders.

This package contains any new features.

- Upgrades the specifications of standard reports in the trial folder.

It upgrades only standard reports that are part of this release.

- Sets permissions on the trial folder and packages for a specific study.

If the file is not run, the packages and the trial folder are only visible to the System Admin user.

Delete obsolete reports

If you are migrating from an InForm 4.6 release, manually delete the following obsolete reports from the study folder:

- CRF Cycle Time Detail by User
- CRF Site Performance Summary by User
- CRF Aging by User
- CRB Completion and Query Summary by User
- Query Volume by Month by User
- Subject Enrollment Status by User

Associate custom reports with the new operational packages

If you are migrating from the InForm 4.6 release:

- 1 Manually associate the custom reports developed in the InForm 4.6 application using the **InForm Trial Management for <trialname>** operational package with the new operational package in InForm 6.0.0.2 **<studyname> InForm Trial Management**.
For each report:
 - a Log in to an InForm study as a user who has the Reporting right and is a member of the Directory Administrators reporting group.
 - b Click **Reports**.
 - c Select **Log On > Log On** to log in to the Cognos 8 Business Intelligence 8.4.1 server as a System Administrator.
 - d In the **User** field, type the user name for the Cognos System Administrator user (for example, crnsysadmin).
 - e In the **Password** field, type the password for the Cognos System Administrator user, and click **OK**.
 - f Select **Actions > Set Properties**.
 - g In the **General** tab, click **Link to a Package**, and select the new **<studyname> InForm Trial Management** package.
 - h Click **OK**.
- 2 Delete the **InForm Trial Management for <studyname>** operational packages from the InForm 4.6 application.

Restart the InForm Model Updater service

- Stop and start the **InForm Model Updater** service to force the immediate clinical model generation.

Revalidate custom reports

- 1 Revalidate the custom reports and, if necessary, modify the reports for the new operational packages.
You revalidate the custom reports by running them and comparing the output to the expected output based on the previous release.
- 2 Enable the `<studyname>` clinical package and study folder.
On the **Properties** page of the clinical package and study folder, deselect **Disable this entry**.
- 3 Stop and start the **InForm Model Updater** service to force the immediate clinical model generation.
- 4 Revalidate all the custom clinical or ad hoc reports for the new clinical model in order to access all the enhancements in this release.
 - If you have access to Report Studio, open the custom clinical report in Report Studio, validate, and save the report.
 - If you do not have access to Report Studio, in the Cognos connection:
 - a Navigate to **Set Properties** for a report, and click **Link to a package**.
 - b Select the **<trial name> clinical package**, and click **OK**.
- 5 Run all the custom clinical reports, and compare the output to the expected output based on the previous release.

Update schedules and jobs

Any user who has scheduled reports must log in to the Reporting and Analysis module and make the necessary changes to the job schedules.

- 1 Select **Launch > Reporting Administration**.
- 2 In the **Status tab**, click **Schedules**.
- 3 Modify and save each listed schedule and job.
 - a Click the drop-down list next to the schedule name.
 - b Click **Modify the Schedule**.
 - c Make any modification (such as changing the start time by one minute), and click **OK**.
- 4 Ensure that study users copy their My Folders content from Public Folder to My Folders.

Delete the study organizational units from the LDAP configuration

Note: This step is not required if you are upgrading from any InForm 6.0 release.

LDAP is not used to authenticate study users in the InForm 6.0.0.2 release.

- Manually delete the study organizational units (OUs) from the Oracle Directory Server using the Oracle Directory Service Control Center.

Running an in-place upgrade of an InForm 6.0 reporting environment

For information on running an in-place upgrade of any InForm 6.0 Reporting environment, see step 10 in *Performing an in-place upgrade for any InForm 6.0 release* (on page 38).

CHAPTER 5

Migrating the InForm Portal software

In this chapter

Prepare the InForm Portal 4.6 database schema for backup.....62

Back up the InForm Portal database and customized files.....63

Migrate the InForm Portal software to a new database instance.....65

Migrate the InForm Portal software in the same database instance.....68

Reset the InForm Portal Home page69

Prepare the InForm Portal 4.6 database schema for backup

If you are migrating an InForm 4.6 study, you must prepare the InForm Portal 4.6 database schema for backup before migrating to the InForm 6.0.0.2 release.

Note: Character set differences between the InForm 4.6 release and the InForm 6.0 releases might cause issues with table column widths. Perform this step to prevent these issues.

- 1 Copy the InForm\InstallSupport folder from the downloaded ISO image to the InForm 4.6 application server.
- 2 Run **Upgrade40toMLInFormPortal.cmd** (located in the InForm\InstallSupport folder) on the InForm 4.6 application server to prepare the InForm 4.6 study schemas for backup, using the following syntax:

```
Upgrade40toMLInFormPortal.cmd <Trial_name> <TNS_Service_Name>
```

where:

- **Trial_Name**—Name of the trial
- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **portaluid**—Owner of the InForm Portal database schema.
- **portalpid**—Password for the owner of the InForm Portal database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

Back up the InForm Portal database and customized files

Back up the InForm Portal 4.6 or 5.0 database and customized files—Oracle 10g

Use this procedure if the InForm Portal software uses an Oracle 10g database.

- 1 Copy the InForm\InstallSupport folder from the downloaded ISO image to the InForm application server.
- 2 Run **ExportDBUpgrade.cmd** (located in the InForm\InstallSupport folder) on the InForm application server to back up the InForm Portal database as a DMP file, using the following syntax:

```
ExportDBUpgrade.cmd <DMP_file> <TNS_Service_Name>
```

where:

- **DMP_file**—Name of the DMP file, using the full path.

Note: Do not include spaces in the DMP file name.

- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **trial_user**—Owner of the InForm Portal database schema.
- **trial_user_pass**—Password for the owner of the InForm Portal database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

Example:

```
ExportDBUpgrade.cmd D:\PortalDB.dmp TRIAL1
```

The ExportDBUpgrade command creates the log file **<DMP_file>.log** in the current directory.

- 3 Back up any InForm Portal files that you customized.

Back up the InForm Portal 4.6, 5.5, or 6.0 database and customized files— Oracle 11g

Use this procedure if the InForm Portal software uses an Oracle 11g database.

- 1 Copy the InForm\InstallSupport folder from the downloaded ISO image to the InForm database server.
- 2 Run **ExportDB.cmd** (located in the InForm\InstallSupport folder) on the InForm database server to back up the InForm Portal database as a DMP file, using the following syntax:

```
ExportDB.cmd <TNS_Service_Name> <Physical_Directory_Path> <DMP_File_Name>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.
- **Physical_Directory_Path**—Location of the directory that you created on the InForm database server for the data pump.

Note: Do not include spaces in the directory path to the DMP file.

- **DMP_File_Name**—Name of the DMP file to export.

Note: Do not include spaces in the DMP file name.

You are prompted for the following values:

- **orasystem_user**—Name of the Oracle database system user.
- **orasystem_user_pass**—Password for the Oracle database system user.
- **trial_user**—Name of the study database schema owner.
- **trial_user_pass**—Password for the owner of the study database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

For example:

```
ExportDB.cmd trial11 /u01/app/oracle/dp_export pfst55.dmp
```

The ExportDB command creates the log file **<DMP_File>.log**.

- 3 Back up any InForm Portal files that you customized.

Migrate the InForm Portal software to a new database instance

- 1 Uninstall the previous version of the InForm Portal software.
 - a Manually remove the InForm Portal DSN entries.
 - b Manually remove the InForm Portal virtual directories.
- 2 Run **DropDB.cmd** (located in the *<InstallationDirectory>\InForm\Bin\DBOra* folder) to create a database user, using the following syntax:

```
DropDB.cmd <TNS_Service_Name>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **pfdbauid**—InForm database admin user.
- **pfdbapid**—InForm database admin password.
- **trial_user**—Owner of the InForm Portal database schema created by the import.
- **trial_user_pass**—Password for the owner of the InForm Portal database schema created by the import.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

Example:

```
DropDB.cmd TRIAL1
```

The DropDB command creates the log file **dropdb.err**.

- 3 Import the InForm Portal 4.6 or 5.0 database schema, and then proceed to step 5.

Note: To import an InForm Portal 5.5 or 6.0 database schema, proceed to step 4.

Run **ImportDBUpgrade.cmd** (located in the *<InstallationDirectory>\InForm\Bin\DBOra* folder) to import the InForm Portal database schema, using the following syntax:

```
ImportDBUpgrade.cmd <DMP_file> <TNS_Service_Name> <From_Portal_UserID>  
<To_Portal_UserID>
```

where:

- **DMP_file**—Name of the DMP file, using the full path.
- **TNS_Service_Name**—TNS name for the database instance.
- **From_Portal_UserID**—Owner of the portal database schema being migrated.
- **To_Portal_UserID**—Owner of the Portal database schema created for the import.

You are prompted for the following values:

- **orasystem_user**—Oracle system user ID.
- **orasystem_user_pass**—Oracle system user password.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is **parameter=value**. There is a new line for each parameter, and there are no spaces on a line.

Example:

```
ImportDBUpgrade.cmd D:\PORTAL46.dmp TRIAL1 PORTAL46UID PORTAL46UID
```

The **ImportDBUpgrade** command creates the log file **<DMP_file>.log**. For the previous example, the log file is **portal46.log**.

- 4 Import the InForm Portal 5.5 or 6.0 database schema.

Run **ImportPortalDB.cmd** (located in the **<InstallationDirectory>\InForm\Bin\DBOra** folder) to import the InForm Portal database schema, using the following syntax:

```
ImportPortalDB.cmd <TNS_Service_Name> <Physical_Directory_Path>  
<DMP_File_Name> <From_Portal_UserID>
```

where:

- **TNS_Service_Name**—TNS name for the database instance.
- **Physical_Directory_Path**—Location of the directory that you created on the InForm database server for the data pump.
- **DMP_File_Name**—Name of the DMP file to import.
- **From_Portal_UserID**—Owner of the portal database schema being migrated.

You are prompted for the following values:

- **orasystem_user**—Oracle database system user.
- **orasystem_user_pass**—Password for the Oracle database system user.
- **trial_user**—Owner of the new portal database schema created by the **ImportPortalDB** command.
- **trial_user_pass**—Password of the owner of the new portal database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is **parameter=value**. There is a new line for each parameter, and there are no spaces on a line.

For example:

```
ImportPortalDB.cmd trial1 /u01/app/oracle/dp_export pfst55.dmp pfst55uid
```

The **ImportPortalDB** command creates the log file **<DMP_file>.log**.

- 5 Run **Upgrade40toMLInFormPortal.cmd** (located in the InForm\InstallSupport folder), using the following syntax:

```
Upgrade40toMLInFormPortal.cmd <Trial_name> <TNS_Service_Name>
```

where:

- **Trial_Name**—Name of the trial or study.
- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **portaluid**—Owner of the InForm Portal database schema.
- **portalpid**—Password for the owner of the InForm Portal database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is parameter=value. There is a new line for each parameter, and there are no spaces on a line.

Example:

```
Upgrade40toMLInFormPortal.cmd PFST55 TRIAL1
```

- 6 Run the Portal Configuration Tool.

Be sure to *deselect* the **Create Portal DB Schema** option, otherwise the Portal Configuration Tool drops the imported database and creates a new one. For more information, see the *Portal Administration Guide*.

- 7 If you are upgrading the InForm Portal from a release earlier than InForm 4.6 SP3 or InForm 5.0 SP2 and you have an existing InForm Portal Home page and images, you must copy them to a new file location after the upgrade is complete.

Move the files from:

- **Homedefault file:** <InstallationDirectory>\InForm\InFormPortal\<locale>\Users (InForm 5.0)
or <InstallationDirectory>\InForm\InFormPortal\Users (InForm 4.6)
- **Images:** <InstallationDirectory>\InForm\InFormPortal\<locale>\Images (InForm 5.0)
or <InstallationDirectory>\InForm\InFormPortal\Images (InForm 4.6)

To:

- **Homedefault file:**
<InstallationDirectory>\InForm\InFormPortal_<studyname>\Views\Display
- **Images:** <InstallationDirectory>\InForm\InFormPortal_<studyname>\Images

- 8 Reset the InForm Portal Home page, if required.

For more information, see *Reset the InForm Portal Home page* (on page 69).

Migrate the InForm Portal software in the same database instance

- 1 Uninstall the previous version of the InForm Portal software.
 - a Manually remove the InForm Portal DSN entries.
 - b Manually remove the InForm Portal virtual directories.
- 2 If you are migrating from any InForm 4.6, InForm 4.7, or InForm 5.0 release, run **Upgrade40toMLInFormPortal.cmd** (located in the *<InstallationDirectory>\InForm\InFormPortal\Scripts\InstallSupport* folder) for each Portal schema, using the following syntax:

Note: This step is not required if you are migrating from any InForm 5.5 or 6.0 release.

```
Upgrade40toMLInFormPortal.cmd <Trial_Name> <TNS_Service_Name>
```

where:

- **Trial_Name**—Name of the study being migrated.
- **TNS_Service_Name**—TNS name for the database instance.

You are prompted for the following values:

- **portaluid**—Owner of the InForm Portal database schema.
- **portalpid**—Password for the owner of the InForm Portal database schema.

You can also pass the parameters by using the **/accountparams:"path_to_parameter_file"** command option.

When specified, this option includes the path to a text file that contains the values required to run the command. The format of the parameter file is **parameter=value**. There is a new line for each parameter, and there are no spaces on a line.

Example:

```
Upgrade40toMLInFormPortal.cmd PFST46 TRIAL1
```

- 3 Run the Portal Configuration Tool.

Be sure to *deselect* the **Create Portal DB Schema** option, otherwise the Portal Configuration Tool drops the imported database and creates a new one. For more information, see the *Portal Administration Guide*.
- 4 Reset the InForm Portal Home page, if required.

For more information, see **Reset the InForm Portal Home page** (on page 69).

Reset the InForm Portal Home page

If the InForm Portal Home page URL has changed as part of the upgrade, reset the home URL for each user.

- 1 Log in to the InForm application as an administrator who can modify user settings.
- 2 Update the Home page for each InForm Portal user to **./portal/default.htm**.

CHAPTER 6

Database and file changes

In this chapter

Database schema changes.....72

Data changes.....79

MedML schema changes81

Base MedML changes85

Design changes.....91

Database and file changes to support the InForm Portal migration.....92

Database schema changes

Changes to the database schema from the InForm 4.6 software to the InForm 4.7 software

The following table lists the database objects that changed from the InForm 4.6 SP2 software to the InForm 4.7 SP0 software.

Database Table	Changes
INF_MSGEDIT_HISTORY	TRANSACTIONCONTENT – LONG to CLOB
PF_COMMENT	COMMENTTEXT – LONG to CLOB
PF_CONTROLDATA	BLOBVALUE – LONG RAW to CLOB. Data is converted from character set AMERICAN_AMERICA.WE8ISO8859P1 to character set AMERICAN_AMERICA.AL32UTF8 to support migration of LATIN 1 characters.
PF_CRFSPATIENTFORM	XMLDATA – LONG RAW to BLOB
PF_EXECUTIONPLAN	EPSCRIPT – LONG to CLOB
PF_EXECUTIONPLANQUEUE	CONTEXTDATA – LONG RAW to BLOB
PF_GROUP	GROUPNAME – Length to 1000
PF_ITEMDATA	ENTEREDREASON – Length from 255 to 1000
PF_LANGUAGE	1 CORE – New Column – Number 2 STUDY – New Column – Number
PF_PATCH	Removed.
PF_PATH	1 PATHSTRING – Length from 255 to 510 2 REFNAME – Length from 255 to 510 3 TYPE – Length from 32 to 64
PF_PATIENT	1 INITIALS – CHAR(3 BYTE) to VARCHAR2(3 CHAR) 2 SCREENINGFAILURE – 255 to 1000 3 ENROLLMENTFAILURE – 255 to 1000 4 ENROLLMENTOVERRIDE – 255 to 1000
PF_PATIENT_CUR	1 INITIALS – CHAR(3 BYTE) to VARCHAR2(3 CHAR) 2 SCREENINGFAILURE – 255 to 1000 3 ENROLLMENTFAILURE – 255 to 1000 4 ENROLLMENTOVERRIDE – 255 to 1000
PF_RESOURCEDATA	RESOURCEDATA – LONG RAW to BLOB

Database Table	Changes
PF_REVISIONHISTORY	REASON – Length from 255 to 500
PF_RIGHTS	RIGHTDESCRIPTION – Length from 63 to 126 bytes
PF_RULES	<ol style="list-style-type: none"> 1 RULEDESCRIPTION – Length from 126 to 512 bytes 2 RULESCRIPT – LONG to CLOB
PF_SITE	<ol style="list-style-type: none"> 1 SITENAME – Length from 255 to 1000 2 SITEMNEMONIC – Length from 32 to 64 3 EMAILADDRESS – Length from 256 to 1000
PF_SITE_CUR	<ol style="list-style-type: none"> 1 SITENAME – Length from 255 to 1000 2 SITEMNEMONIC – Length from 32 to 64 3 EMAILADDRESS – Length from 256 to 1000
PF_SYNCLOG	SYNCLOG – LONG RAW to BLOB
PF_USER	<ol style="list-style-type: none"> 1 FIRSTNAME – Length to 510 2 LASTNAME – Length to 510 3 EMAILADDRESS – Length to 1000 4 COUNTRY – New Column – Varchar2(510)
PF_USERPROFILE	<ol style="list-style-type: none"> 1 LANGUAGEID – New Column Number 2 DISPLAYNAME – New Column Varchar2(255)

Changes to the database schema from the InForm 4.7 software to the InForm 5.0 software

The following table lists the database objects that changed from the InForm 4.7 SP0 software to the InForm 5.0 SP0 software.

Database Table	Change
PF_TextControl	Added CHARSET - New Column Number

Changes to the database schema from the InForm 5.0 software to the InForm 5.5 software

The following table lists the database objects that changed from the InForm 5.0 SP0 software to the InForm 5.5 SP0 software.

Database Table	Column	Column Type
PF_REVISIONHISTORY	REVIEWSCHEMAUPDATED – New Date Column	
PF_SVSUBJECTPOOL (new table)	SVSUBJECTPOOLID	NUMBER NOT NULL
	SVSUBJECTPOOLREVISIONNUMBER	FLOAT NOT NULL
	PATIENTID	NUMBER NOT NULL
	SITEID	NUMBER NOT NULL
	POOLNUM	NUMBER NOT NULL
PF_SVSUBJECTPOOL_CUR (new table)	SVSUBJECTPOOLID	NUMBER NOT NULL
	SVSUBJECTPOOLREVISIONNUMBER	FLOAT NOT NULL
	PATIENTID	NUMBER NOT NULL
	SITEID	NUMBER NOT NULL
	POOLNUM	NUMBER NOT NULL
PF_SVCriticalForms (new table)	SVCriticalFORMSID	NUMBER NOT NULL
	SVCriticalFORMSREVISIONNUMBER	FLOAT NOT NULL
	PAGEID	NUMBER NOT NULL
	SITEID	NUMBER NOT NULL
	ISCRITICAL	NUMBER NOT NULL
PF_SVCriticalForms_CUR (new table)	SVCriticalFORMSID	NUMBER NOT NULL
	SVCriticalFORMSREVISIONNUMBER	FLOAT NOT NULL
	PAGEID	NUMBER NOT NULL
	SITEID	NUMBER NOT NULL
	ISCRITICAL	NUMBER NOT NULL

Changes to the database schema from the InForm 5.5 software to the InForm 6.0 software

The following table lists the database objects that changed from the InForm 5.5 SP0 software to the InForm 6.0 software.

Database Table	Column	Column Type
PF_COMMENT	REVIEWER	New Number Column
	REVIEWDATE	New Date Column
PF_CONTROL	ISCOLLAPSIBLE	New Number Column
PF_ITEM	SDVCRITICAL	New Number Column
PF_ITEMCONTEXT	PREDEFINED	New Number Column
	BLANK	New Number Column
PF_USERPROFILE	STRVALUE	Varchar2 to CLOB
PF_PATH	FORMVERSION	New Number Not Null
PF_ITEMSET(new table)	ITEMSETID	Number Not Null
	ITEMSETREVISIONNUMBER	FLOAT Not Null
	ITEMSETREFNAME	VARCHAR2(4000) Not Null
	INITIALROWCOUNT	NUMBER Not Null
	ROWCOUNTFIXED	NUMBER Not Null
PF_ITEMSETCELL(new table)	ITEMSETID	Number Not Null
	ITEMSETREVISIONNUMBER	Float Not Null
	ITEMID	Number Not Null
	ITEMSETROW	Number Not Null
	PREDEFINEDELEMENTID	Number Not Null
	ISBLANK	Number Not Null
PF_SVCRITICALITEMS(new table)	SVCRITICALITEMSID	Number Not Null
	SVCRTITEMSREVISIONNUMBER	Float Not Null
	SITEID	Number Not Null
	PAGEID	Number Not Null
	ITEMID	Number
	ISCRITICAL	

Database Table	Column	Column Type
PF_SVCRITICALITEMS_CUR(Not Null)	SVCRITICALITEMSID	Number Not Null
	SVCRTTITEMSREVISIONNUMB ER	Float Not Null Number Not Null
	SITEID	Number Not Null
	PAGEID	Number Not Null
	ITEMID	Number
	ISCRITICAL	
PF_XMLSUBMITDOCUMENT STATUS(new table)	XMLSUBMITDOCUMENTID	VARCHAR(255)
	XMLSUBMITDOCUMENTTYPE	Number Not Null
	RECIEVEDATETIME	TIMESTAMP Not Null
	PROCESSINGSTARTDATETIME	TIMESTAMP Not Null
	PROCESSINGSTOPDATETIME	TIMESTAMP
	TOTALSUBMITS	Number Not Null
	TOTALTRANSACTIONS	Number Not Null
	RULESEENABLED	Number Not Null
	STOPONERROR	Number Not Null
	SUBMITTOLOCKEDORFROZE NFORMS	Number Not Null
	USERNAME	VARCHAR(255) Not Null
	TOTALSUBMITSSUCCEDED	Number
	TOTALSUBMITSFAILED	Number
	TOTALSUBMITSNOTEXECUTE D	Number Number
	COMMITTEDTRANSACTIONS	Number

Database Table	Column	Column Type
PF_XMLSUBMITCOMMAND DATA (new table)	XMLSUBMITDOCUMENTID	VARCHAR(255) NOT NULL
	INFORMCOMMAND	NUMBER NOT NULL
	SUBJECTIDENTIFIER	VARCHAR(255) NOT NULL
	LOCATIONOID	VARCHAR(255) NOT NULL
	METADATAVERSIONOID	VARCHAR(255) NOT NULL
	STUDYEVENTOID	VARCHAR(255)
	STUDYEVENTREPEATKEY	VARCHAR(255)
	FORMOID	VARCHAR(255)
	FORMREPEATKEY	VARCHAR(255)
	SUBMITINDEX	VARCHAR(255)
	TRANSACTIONINDEX	VARCHAR(255)
	LINENUMBER	NUMBER NOT NULL
	LINEPOSITION	NUMBER NOT NULL
	SUBMITSTARTDATETIME	NUMBER NOT NULL
	SUBMITSTOPDATETIME	NUMBER NOT NULL
	ACTIVITYGUID	TIMESTAMP NOT NULL
	ERRORMESSAGE	TIMESTAMP
	RESULTCODE	VARCHAR(255)
PF_XMLSUBMITDOCUMENT DATA(new table)	XMLSUBMITDOCUMENTID	VARCHAR(255) NOT NULL
	XMLSUBMITDOCUMENT	CLOB NOT NULL
	XMLSUBMITDOCUMENTFILE NAME	VARCHAR(255)
PF_BATCHRUNRULES (new table)	BATCHRUNRULEID	NUMBER NOT NULL
	BATCHRUNRULEREVISIONNU MBER	FLOAT NOT NULL
	PFTimeStamp	DATE NULL
	STATUS	NUMBER NOT NULL
	DETAILS	VARCHAR2(255 CHAR)

Database Table	Column	Column Type
PF_BATCHRUNRULEDETAILS (new table)	BATCHRUNRULEID	NUMBER NOT NULL
	PATIENTID	NUMBER NOT NULL
	RULEID	NUMBER NOT NULL
	RULECONTEXTID	NUMBER NOT NULL
	PFTimeStamp	DATE NULL
	STATUS	NUMBER
	DETAILS	VARCHAR2(255 CHAR)
PF_PATIENTTRANSFERSTATUS (new table)	TRANSACTIONID	NUMBER NOT NULL
	ARCHIVEREVISIONNUMBER	FLOAT NOT NULL
	PATIENTID	NUMBER NOT NULL
	PATIENTREVISIONNUMBER	FLOAT NOT NULL
	FROMSITEID	NUMBER NOT NULL
	TOSITEID	NUMBER NOT NULL
	USERID	NUMBER NOT NULL
	NEWPATIENTNUMBER	VARCHAR(255)
	REASON	VARCHAR(255)
	PFTimeStamp	DATE
	STATUS	NUMBER NOT NULL
	DETAILS	VARCHAR(255)
PF_ITEMSETROWCONTEXT (new table)	CONTEXTID	NUMBER NOT NULL
	ITEMSETREPEATKEY	VARCHAR2(255 CHAR) Not Null

Data changes

Data changes from the InForm 4.6 software to the InForm 4.7 software

The following table lists the data changes from the InForm 4.6 SP2 software to the InForm 4.7 SP0 software.

Table	Changes
PF_INFORM	Update DBVersion to InForm 4.7 .
PF_DBUID	Update GUID to English (United States) where DBUID is 2.
PF_LANGUAGE	<ol style="list-style-type: none"> 1 Change LanguageName to en-US from English. 2 Set CORE to 1 where Language is en-US. 3 Set STUDY to 0 where Language is en-US. 4 Add Japanese language ja-JP.
PF_DBUID	Add Japanese Language ID.
PF_ENTITYMAPPING	Add Japanese Language LCID.
PF_USERPROFILE	Add DisplayName entries Not Complete , Not Verified , and Ready to Verify for English and Japanese languages.

Data changes from the InForm 4.7 software to the InForm 5.0 software

None.

Data changes from the InForm 5.0 software to the InForm 5.5 software

The following table lists the data changes from the InForm 5.0 SP0 software to the InForm 5.5 SP0 software.

Table	Changes
PF_INFORM	Update DBVersion to InForm 5.5 .
PF_TABLES	The following new values were added: 1 PF_CHAPTER_CUR 2 PF_IM_CHAPTER 3 PF_IM_CHAPTER_CUR 4 PF_IM_PAGE 5 PF_PATIENTNUMEBR_CUR 6 PF_PATIENT_CUR 7 PF_QUERY_LISTING 8 PF_RESOURCEDATA_CUR 9 PF_SITEVOLUMEEDITION_CUR 10 PF_SITE_CUR 11 PF_SUBJECTVECHAPTER_CUR 12 PF_USERPROFILE

Data changes from the InForm 5.0 software to the InForm 5.5 software

The following table lists the data changes from the InForm 5.5 SP0 software to the InForm 6.0 software.

Table	Changes
PF_INFORM	Update DBVersion to InForm 6.0
PF_TABLES	The following new values were added: 1 PF_CHAPTER_CUR 2 PF_IM_CHAPTER

MedML schema changes

MedML schema changes from the InForm 4.6 software to the InForm 4.7 software

The following table lists the MedML database schema objects that changed from the InForm 4.6 SP2 software to the InForm 4.7 SP0 software.

Category	Changes
New attributes	1 TypeNonEmptyString
	2 TypeTranslatableAttribute
	3 TypeDATASOURCE
New elements	1 TRANSLATION
	2 TRANSLATIONS
Changed attribute	TypeLanguage
Elements changed to have optional TRANSLATIONS element	1 DOCUMENTATION
	2 DOCBODY
	3 CALCLATEDCONTROL
	4 ITEM
	5 SECTION
	6 UNIT
	7 PFELEMENT
	8 SIMPLECONTROL
	9 TEXTCONTROL
	10 RADIOCONTROL
	11 DATETIMECONTROL
	12 FORM
	13 PULLDOWNCONTROL
	14 CHECKBOXCONTROL
	15 LISTBOXCONTROL
	16 GROUPCONTROL
	17 FORMSET
	18 SIGNATUREGROUP

Category	Changes
Changed elements	1 USER: New required attributes PRODUCTLOCALE and STUDYLOCALE.
	2 SYSCONFIG: TYPE is now an optional attribute.
	3 SITE: New optional attributes STUDYLOCALE and USERNAMEORDER.
	4 SIGNATUREGROUP: New optional attribute UPDATE.

MedML schema changes from the InForm 4.7 software to the InForm 5.0 software

The following table lists the MedML database schema objects that changed from the InForm 4.7 SP0 software to the InForm 5.0 SP0 software.

Element	Changes
TEXTCONTROL MedML	Added CHARSET

MedML schema changes from the InForm 5.0 software to the InForm 5.5 software

The following table lists the MedML database schema objects that changed from the InForm 5.0 SP0 software to the InForm 5.5 SP0 software.

Category	Changes
Changed element: SITE	1 New optional attribute UPDATE to support updating existing site objects.
	2 New optional attribute SVAUTOSELECTRATE to support the Partial SV feature.
	3 New optional SVFIRSTNSUBJECTS to support Partial SV.
	4 Attribute NAME changed to REQUIRED from OPTIONAL.

MedML schema changes from the InForm 5.5 software to the InForm 6.0 software

The following table lists the MedML database schema objects that changed from the InForm 5.5 SP0 software to the InForm 6.0 software.

Category	Changes
Changed element: ITEM	New optional attribute SDVCRITICAL of type BOOL
Changed element: TEXTCONTROL	New optional attribute ISCOLLAPSIBLE of type BOOL
Changed element: RADIOCONTROL	New optional attribute ISCOLLAPSIBLE of type BOOL
Changed element: DATETIMECONTROL	New optional attribute ISCOLLAPSIBLE of type BOOL
Changed element: PULLDOWNCONTROL	New optional attribute ISCOLLAPSIBLE of type BOOL
Change element: CHECKBOXCONTROL	New optional attribute ISCOLLAPSIBLE of type BOOL
Changed element: GROUPCONTROL	New optional attribute ISCOLLAPSIBLE of type BOOL
Changed element: USER	<ol style="list-style-type: none"> 1. New optional attribute SALT of type string 2. New optional attribute HASHTYPE of type string 3. New optional attribute TRUSTPEUID of type string
New element: USERNAMEUPDATE	<ol style="list-style-type: none"> 1. Required attribute CURRENTUSERNAME of type TypeUSERNAME 2. Required attribute NEWUSERNAME of type USERNAME 3. Optional attribute TRUSTPEUID of type USERNAME
New element: ITEMSETCELL	<ol style="list-style-type: none"> 1. Required attribute ITEMREFNAME of REFNAME 2. Required attribute ITEMSETROW of type positiveInteger 3. Optional attribute PREFEDEFINEDVALUEREFFNAME of type REFNAME 4. Optional attribute ISBLANK of type BOOL
Changed element: ITEMSET	<ol style="list-style-type: none"> 1. Element ITEMSETCELL added to the sequence 2. New optional attribute SDVCRITICAL of type BOOL 3. New optional attribute UUID of type GUID 4. New optional attribute INITIALROWCOUNT of type nonNegativeInteger 5. New optional attribute ROWCOUNTFIXED of type BOOL

Category	Changes
New element: SVCRITICALFORM	<ol style="list-style-type: none">1. Required attribute FORMREFNAME of type REFNAME2. Optional attribute CRITICAL of type BOOL3. Optional attribute SITEGUID of type GUID4. Optional attribute ITEMNEMONIC of type MNEMONIC5. Optional attribute GUID of type GUID6. Optional attribute REVISION of type REVISION
New element SVCRITICALITEM	<ol style="list-style-type: none">1. Required attribute CRITICAL of type string2. Required attribute ITEMREFNAME of type REFNAME3. Optional attribute GUID of type GUID4. Optional attribute REVISION of type REVISION5. Optional attribute FORMREFNAME of type REFNAME6. Optional attribute SITEGUIDREF of type GUIDREF

Base MedML changes

Base MedML changes from the InForm 4.6 software to the InForm 4.7 software

The following table lists the Base MedML that changed from the InForm 4.6 SP2 software to the InForm 4.7 SP0 software. All files are located in the folder <InstallationDirectory>\InForm\XMLBase.

File	Description	Change
Core_Date.xml	Core Date Formats.	1 Added visit formats for Visit Calc, Year suffix, and Month suffix. 2 Added Date Formats for Japanese.
Core_doc.xml	DOCUMENTATION objects.	Added translations for DOCNAME and HELPTTEXT attributes.
Core_grouptypes.xml	THINGTYPE properties.	Added properties for SiteStudyLocale and SiteUserNameOrder.
Core_Patch.xml	Patch things.	Removed from XMLBase
Core_PatchAdmin.xml	Patch Admin settings.	Removed from XMLBase
Core_Randomization.xml	Randomization section.	Added translation to objects.
Core_RightsGroups.xml	RightsGroups.	Removed rights Manage Software Patch and View Software Patch Status rights from SysAdmin Rights Group RightsGroup.
Core_systemConfig.xml	System config data.	Removed value for NumCharsCRFLength . Added value for USERNAMEORDER . Note: Installing this file into an existing study changes the existing system configuration settings.
Core_SystemConfigThing.xml	THINGTYPE properties.	Added properties UserProductLocale and UserStudyLocale .
Core_ValuesMetaData.xml	DOV section.	Added translation to objects.
Sysform_AdminDBSvr.xml	AdminDB server context form.	Added translation to objects.
Sysform_AdminDBTrl.xml	AdminDB study context form.	Added translation to objects.
Sysform_AdminUserNew.xml	User Admin form.	Added translation to objects.

File	Description	Change
Sysform_AdvancedConnectionSettings.xml	Advanced Connection Settings.	Added translation to objects.
Sysform_AssignUserToGroup.xml	User group edit form.	Added translation to objects.
Sysform_ClearCRF.xml	Clear CRF form.	Added translation to objects.
Sysform_Comment.xml	Comment form.	Added translation to objects.
Sysform_Conflict.xml	Conflict resolution.	Added translation to objects.
Sysform_ConnectionSettings.xml	Connection settings.	Added translation to objects.
Sysform_DefaultAdvancedConnectionSettings.xml	Default Advanced settings.	Added translation to objects.
Sysform_DefaultConnectionSettings.xml	Sync Default settings.	Added translation to objects.
Sysform_DictType.xml	Dictionary properties form.	Added translation to objects.
Sysform_EditItem.xml	Edit Item form.	Added translation to objects.
Sysform_EventAdd.xml	Add event form.	Added translation to objects.
Sysform_FAQ.xml	FAQ form.	Added translation to objects.
Sysform_Group.xml	Add group form.	Added translation to objects. The file is obsolete.
Sysform_GroupAddProp.xml	Add group properties.	Added translation to objects.
Sysform_GroupEditProp.xml	Edit group properties.	Added translation to objects.
Sysform_GroupSelect.xml	User group association.	Added translation to objects.
Sysform_GroupSelect.xml	User group association.	Added translation to objects.
Sysform_GroupType.xml	Reporting group form.	Added translation to objects.
Sysform_ItemEditReasons.xml	Edit item reason.	Added translation to objects.
Sysform_Listings.xml	Internal listings form.	Added translation to objects.
Sysform_Password.xml	Password entry form.	Added translation to objects.
Sysform_Patch.xml	Patch form.	Removed.
Sysform_PatchAdminThings.xml	Patch form.	Removed.
Sysform_PatchForms.xml	Patch form.	Removed.
Sysform_PatientSiteChange.xml	Patient site change form.	Added translation to objects.
Sysform_Query.xml	Query form.	Added translation to objects.
Sysform_Query2.xml	Query form.	Added translation to objects.

File	Description	Change
Sysform_RightsEdit.xml	Right edit form.	Added translation to objects.
Sysform_Rule.xml	Rule edit form.	Added translation to objects.
Sysform_RuleItemEdit.xml	Rule item edit form.	Added translation to objects.
Sysform_Site.xml	Site edit form.	Added translation to objects.
Sysform_synchronization.xml	Sync file transfer.	Added translation to objects.
Sysform_SystemConfig.xml	System Configuration.	Added translation to objects.
Sysform_SystemInfo.xml	System Info.	Added translation to objects.
Sysform_UserPswrd.xml	User Password form.	Added translation to objects.
Sysform_UserType.xml	User Information form.	Added translation to objects.
sysVolume.xml	Version Information.	Changed version.
sysVolumeAdmin.xml	System forms.	Added translation to objects.
QueryFormText.xml	Query form.	New file.

Base MedML changes from the InForm 4.7 software to the InForm 5.0 software

The following table lists the Base MedML that changed from the InForm 4.7 SP0 software to the InForm 5.0 SP0 software. All files are located in the folder <InstallationDirectory>\InForm\XMLBase.

File	Description	Change
sysform_SystemConfig.xml	System Configuration Form	<ol style="list-style-type: none"> 1 Added a read-only display indicating if forms display traffic light icons. 2 Increased the E-mail address string max size to 255 characters.
core_SystemConfig.xml	Default values for system configuration parameters.	<ol style="list-style-type: none"> 1 Changed QueryMaxLength from 80 to 350 characters. 2 Added DISPLAYTRAFFICLIGHTS and set the default to false (0). <p>Note: Installing this file into an existing study changes the existing system configuration settings.</p>

File	Description	Change
core_SystemConfigThing.xml	Thing types, part of Thing World: The InForm application's means of categorizing and storing property values.	Added the DISPLAYTRAFFICLIGHTS thing type.
core_GroupTypeStrings.xml	Translations for thing type descriptions.	New file.
sysVolumeAdmin.xml	System Forms	Added LANGUAGE="en-US"
sysform_SystemInfo.xml	Copyright info	Change year from 2008 to 2009 in copy rights
sysform_Site.xml	Site details page	Updated patient to subject .
sysform_RightsEdit.xml	Rights Group details page	Updated patient to subject .
sysform_Password.xml	Update user password page.	Updated translations.
sysform_PatientSiteChange.xml	Patient Site Change page.	Patient changed to Subject .
core_Date.xml	Date and Time formats	Date and Time separators
XmlBaseEnd.xml	End Of Base XML	New file.

Base MedML changes from the InForm 5.0 software to the InForm 5.5 software

The following table lists the Base MedML that changed from the InForm 5.0 SP0 software to the InForm 5.5 SP0 software. All files are located in the folder *<InstallationDirectory>\InForm\XMLBase*.

File	Description	Changes
Core_Rights.xml	Rights	<p>Added the Assign Review States and View Review States rights, which provide access to the custom review states in the Data Viewer.</p> <p>Added the Change SV Site Settings right, which allows access to the Partial SV configuration options.</p>
Core_RightsGroups.xml	Rights Groups	Changed rights groups System Creator Group and InForm Server Group .
Core_SystemConfig.xml	System config values	<ol style="list-style-type: none"> Added value for Authentication Component. Changed value for UNIQUEPATIDSWTCH. <p>Note: Installing this file into an existing study changes the existing system configuration settings.</p>

File	Description	Changes
Core_SystemConfigThin.g.xml	System Config things	Added AUTHENTICATIONCOMPONENT.
Sysform_RightsEdit.xml	Rights Edit form	Added Data Viewer section. Note: During migration to the InForm 5.5 release, the InForm 5.0 user right Data Export Listings is changed to Access Data Viewer / Data Export Listings .
Sysform_SystemConfig.xml	System Config form	Added Authentication Component.
sysVolume.xml	System form version	Change in version.
Reports/Report.xml	Reports xml	Removed.

Base MedML changes from the InForm 5.5 software to the InForm 6.0 software

The following table lists the Base MedML that changed from the InForm 5.5 SP0 software to the InForm 6.0 software. All files are located in the folder *<InstallationDirectory>\InForm\XMLBase*.

File	Description	Changes
Clinical_Users.xml	Users	Ctcoding and ctvalidation users inactivated and terminated.
Core_MercUsers.xml	Users	File removed
Core_ReporingUser.xml	Rights Group	Added "Reports Only Rights Group"
Default_Users.xml	Users	File removed
Core_MercRightGroups	RightsGroup	File removed
Core_users.xml	Users	Removed UM,UA,InFormUser,pfarchuser,sysadmin users.
Core_usertype.xml	Users	Added property "UserTrustPEUID"
Sysvolume.xmltmpl	Study Version	Changed InForm version format Removed sysform_FAQ.xml
Sysform_UserType.xml	User form	Added UserSalt,UserHashType and UserTrustPEUID items
upgradeAdminConfigAttributes.xml	System Config	Removed Authentication Component
Default_rightsgroups.xml	Rights Groups	Removed deleted users from default rights groups membership

File	Description	Changes
Core_Rights.xml	Rights	Added new rights
Core_RightsGroups.xml	Rights Groups	Removed deleted uses from rights groups membership
Core_systemconfig.xml	System Config	Removed authentication component
core_SystemConfigThing.xml	System Config	Removed authentication component Added TrialOrgID and PfReportingUserName properties
UpgradeTrialResources.rsp	RSP file	Removed reference to files sysform_FAQ.xml, core_MercUsers.xml,, core_MercRightsGroups.xml
sysform_RightsEdit.xml	Rights	Added rights
Sysform_systemconfig.xml	System Config	Removed authentication component Added Trial OrgID and Reporting user name
sysform_Password.xml	Password Recovery	Automatic password recovery related changes.
sysform_ItemsetRowComment.xml	Custom form for itemset row comment	New file

Design changes

Design changes for the InForm 6.0 software

The following table lists the files that were created or changed in the InForm 6.0 software to support the migration procedure.

Path	Description	Change
<InstallationDirectory>\InForm\Bin\DBOra\UpgradeTrialSchema.cmd	Command file for database upgrade.	Added command to run the sql file schemaUpgrade_5.6.sql.
<InstallationDirectory>\InForm\Bin\DBOra\UpgradeTrialSchema.sql	SQL file for database upgrade.	Added command to support schema upgrade of 5.5 studies.
<InstallationDirectory>\InForm\Bin\DBOra\SchemaUpgrade_5.6.sql	New SQL file for upgrading schema from the InForm 5.5 software to the InForm 6.0 software.	File contains schema changes from the InForm 5.5 software to the InForm 6.0 software..
<InstallationDirectory>\InForm\Bin\DBOra\TestInFormVersion.sql	SQL file to verify the InForm software version	Added code to verify the InForm 5.5 version.
<InstallationDirectory>\InForm\Bin\DBOra\utilityprocedures.sql	SQL file that contains procedures to support InForm migration.	Added new procedures like XmlSubmitPurgeRecords and ConvertNullableStringToClob .
<InstallationDirectory>\InForm\Bin\DBOra\UpgradeResources.cmd	Command file for upgrading resources.	Added changes to support 5.5 to 6.0 upgrade
<InstallationDirectory>\InForm\XMLBase\UpgradeTrialResources.rsp	Contains resource and BaseXML files to be installed.	Added sysform_ItemsetRowComment.xml and removed upgradeAdminConfigAttributes.xml.
<InstallationDirectory>\InForm\Bin\DBOra\clearcontrolpaths.sql	Sql file to delete PF_PATH records	Script to delete PF_PATH records
<InstallationDirectory>\InForm\Bin\DBOra\ImportPortaldb.cmd	Command to import 5.5 Portal schema	Portal import using data pump
<InstallationDirectory>\InForm\Bin\DBOra\TerminateUsers.cmd	Command to terminate system users and blank out invalid email addresses	Contains commands to generate the required MedMLs and cook them

Database and file changes to support the InForm Portal migration

Database schema changes for the InForm Portal software

Database Table	Changes
PF_CATEGORIES	LANGUAGEID – New Column Varchar2(50)
PF_SUBCAT	LANGUAGEID – New Column Varchar2(50)
PF_DOC_REF	REFNAME – Column changed to Varchar2(150) from Varchar2(50)
PF_CATEGORIES	MAINURL- Column changed to Varchar2(255) from Varchar2(250)
P_CONFIG	DSN – Column changed to Varchar2(255) from Varchar2(75)

Data changes for the InForm Portal

Database Table	Changes
PF_CATEGORIES	Update all rows to have LanguageID value to en-US.
PF_SUBCAT	Update all rows to have LanguageID value to en-US.
PF_CONFIG	DSNUID changed to null
PF_CONFIG	DSNPWD changed to null

Additions to the PF_CATEGORIES table for the InForm Portal software

The following table lists the records that were added to the PF_CATEGORIES table to support the Japanese language in the InForm 6.0 release. The ID value is calculated based on existing records in the database and inserted with each record.

CATEGORYNAME	Description	ACTIVE	TAB ORDER	MAIN URL	STR 1	SECTION	SPONSOR	LANGUAGE ID
フォーラム	トライアルに関する様々なトピックについてのディスカッション。アイデアや意見を同僚と共有することができます。	0	0	/webx_trial Name /	50	0	0	ja-JP

CATEGORYNAME	Description	ACTIVE	TABORDER	MAIN URL	STR 1	SECTION	SPONSOR	LANGUAGE ID
ヘルプ	Portalやスタディに関連するヘルプ	0	0	help.asp	42	0	0	ja-JP
タブ1	トライアル関連資料	1	3	view.asp		1	0	ja-JP
ホーム	InForm Portal ホーム	1	2	main.asp		0	0	ja-JP
管理	管理ポータル	1	99	portal/editnav.asp		0	0	ja-JP
コンタクト	主な問合せ先一覧	0	0	contact.asp		0	0	ja-JP
管理者設定	システム設定	1	1	editnav.asp		0	0	ja-JP

Records inserted into the PF_SUBCAT table for the InForm Portal software

The following table lists the records that were inserted into the PF_SUBCAT table for the InForm 6.0 software for the InForm Portal software. The ID value is calculated based on existing records in the database and inserted with each record.

CATEGORYID	SUBNAME	DESCRIPTION	ACTIVE	TABORDER	LANGUAGEID
PF_CATEGORY ID for help.asp	FAQ	"	1	0	ja-JP
PF_CATEGORY ID for main.asp	ホーム	"	1	0	ja-JP

File changes to support the InForm Portal software migration

The following table lists the files changed to support the InForm Portal migration to the InForm 6.0 software.

Path	Description	Changes
<InstallationDirectory>\InForm\InFormPortal\Scripts\InstallSupport\Upgrade4toMLInFormPortal.cmd	Command file for InForm Portal upgrade.	Command to call upgrade4toMLInFormPortal.sql.
<InstallationDirectory>\InForm\InFormPortal\Scripts\InstallSupport\Upgrade\Upgrade4toMLInFormPortal.sql	SQL file for InForm Portal upgrade.	SQL scripts include script to add LANGUAGEID columns to PF_CATEGORIES and PF_SUBCAT tables, and add default data to support Japanese language.